



US005795725A

United States Patent [19][11] **Patent Number:** **5,795,725****Buechler et al.**[45] **Date of Patent:** **Aug. 18, 1998**[54] **METHODS FOR THE ASSAY OF TROPONIN I AND T AND SELECTION OF ANTIBODIES FOR USE IN IMMUNOASSAYS**[75] **Inventors:** **Kenneth F. Buechler; Paul H. McPherson**, both of San Diego, Calif.[73] **Assignee:** **Biosite Diagnostics Incorporated**, San Diego, Calif.[21] **Appl. No.:** **423,582**[22] **Filed:** **Apr. 18, 1995**[51] **Int. Cl.⁶** **G01N 33/553**[52] **U.S. Cl.** **435/7.1; 435/7.2; 435/7.5; 435/7.71; 435/7.72; 435/7.8; 435/7.9; 435/7.92; 435/7.93; 435/7.94; 435/7.95; 435/960; 435/962; 435/971; 435/973; 436/518; 436/523; 436/528; 436/535; 436/538; 436/811; 436/819**[58] **Field of Search** **435/7.1, 7.2, 7.5, 435/7.72, 7.8, 7.9, 7.92, 7.93, 7.94, 7.95, 960, 962, 971, 973, 7.71; 436/518, 523, 528, 535, 538, 811, 819**[56] **References Cited****U.S. PATENT DOCUMENTS**

4,933,275	6/1990	Wands et al.	435/7
5,051,356	9/1991	Warren, III et al.	435/7.54
5,143,852	9/1992	Valkirs et al.	436/501
5,210,020	5/1993	Kondo et al.	435/7.94
5,246,001	9/1993	Ohman et al.	128/630
5,290,678	3/1994	Jackowski	435/7.4
5,447,846	9/1995	Shinoki et al.	435/7.92
5,604,105	2/1997	Jackowski	435/7.4

FOREIGN PATENT DOCUMENTS

2275774	9/1994	United Kingdom
9427156	11/1994	WIPO

OTHER PUBLICATIONSAnderson et al., "Troponin T Isoform Expression in Humans." *Circulation Research* 69(4):1226-1233 (1991).Ball et al., "Isoform Specific Interactions of Troponin I and Troponin C Determine pH Sensitivity of Myofibrillar Ca²⁺ Activation." *Biochemistry* 33:8464-8471 (1994).Bhayana et al., "Discordance Between Results for Serum Troponin T and Troponin I in Renal Disease." *Clin. Chem.* 41(2):312-317 (1995).Bhavsar et al., "Developmental expression of troponin I isoforms in fetal human heart." *FEBS Letters* 292:(1,2):5-8 (1991).Blechner et al., "4Ca²⁺-Troponin C Forms Dimers in Solution at Neutral pH that Dissociate upon Binding Various Peptides: Small-Angle X-ray Scattering Studies of Peptide-Induced Structural Changes." *Biochemistry* 31:11326-11334 (1992).Borrebaeck, *Antibody Engineering: A Practical Guide*, W.H. Freeman and Company, New York (1992).Burtnick et al., "The isolation and characterization of the tropomyosin binding component (TN-T) of bovine cardiac troponin." *Canadian Journal of Biochemistry* 54(5):546-553 (1976).Byers and Kay, "Hydrodynamic Properties of Bovine Cardiac Troponin-I and Troponin-T." *J. Biol. Chem.* 258:2951-2954 (1983).Collins et al., "Early Diagnosis of Acute Myocardial Infarction with Use of a Rapid Immunochemical Assay of Creatine Kinase MB Isoenzyme." *Clin. Chem.* 39:1725-1728 (1993).Adams et al., "Comparable Detection of Acute Myocardial Infarction by Creatine Kinase MB Isoenzyme and Cardiac Troponin I." *Clin. Chem.* 40(7):1291-1295 (1994).Cummins et al., "Cardiac-specific troponin-I radioimmunoassay in the diagnosis of acute myocardial infarction." *American Heart Journal* 113(6):1333-1344 (1987).Cummins and Perry, "Troponin I from Human Skeletal and Cardiac Muscles." *Biochem. J.* 171:251-259 (1978).Harlow and Lane, *Antibodies: A Laboratory Manual*, Cold Spring Harbor Laboratory (1988).Horwitz et al., "Interaction of Troponin Subunits: The Interaction Between the Inhibitory and Tropomyosin-Binding Subunits." *J. Biol. Chem.* 254(2):350-355 (1979).Huse et al., "Application of a Filamentous Phage pVIII Fusion Protein System Suitable for Efficient Production, Screening, and Mutagenesis of F(ab) Antibody Fragments." *J. Immunology* 149:3915-3920 (1992).Katus et al., "Intracellular Compartmentation of Cardiac Troponin T and its Release Kinetics in Patients with Reperfused and Nonreperfused Myocardial Infarction." *American Journal of Cardiology* 67:1360-1367 (1991).Katus et al., "Diagnostic Efficiency of Troponin T Measurements in Acute Myocardial Infarction." *Circulation* 83(3):902-912 (1991).Leavis and Kraft, "Calcium Binding to Cardiac Troponin C^{1,2}." *Archives of Biochemistry and Biophysics* 186(2):411-415 (1978).Liao et al., "Coupling of Calcium to the Interaction of Troponin I with Troponin C from Cardiac Muscle." *Biochemistry* 33:12729-12734 (1994).Liao et al., "Time-resolved tryptophan emission study of cardiac troponin I." *Biophys. J.* 63:986-995 (1992).MacGeoch et al., "The human cardiac troponin I locus: assignment to chromosome 19p13.2-19q13.2." *Human Genet.* 88:101-104 (1991).Metzger et al., "Skeletal troponin C reduces contractile sensitivity to acidosis in cardiac myocytes from transgenic mice." *Proc. Natl. Acad. Sci. USA* 90:9036-9040 (1993).

(List continued on next page.)

Primary Examiner—Christopher L. Chin
Assistant Examiner—Bao-Thuy L. Nguyen
Attorney, Agent, or Firm—Lyon & Lyon LLP

[57] **ABSTRACT**

Assays and antibodies are disclosed for the detection and quantitation of cardiac specific troponin I and troponin T in body fluids as an indicator of myocardial infarction. Since troponin I and T exist in various conformations in the blood, the ratios of the monomeric troponin I and T and the binary and ternary complexes may be related to the metabolic state of the heart. More specifically, immunoassays for determining the presence or amount of free, binary and ternary complexes of troponin I and T are claimed.

18 Claims, 11 Drawing Sheets

OTHER PUBLICATIONS

- Ngai and Hodges. "Biologically Important Interactions between Synthetic Peptides of the N-terminal Region of Troponin I and Troponin C." *J. Biol. Chem.* 267:15715-15720 (1992).
- Olah and Trewhella. "A Model Structure of the Muscle Protein Complex $4Ca^{2+}$ -Troponin C-Troponin I Derived from Small-Angle Scattering Data: Implications for Regulation." *Biochemistry* 33:12800-12806 (1994).
- Olah et al.. "Troponin I Encompasses an Extended Troponin C in the Ca^{2+} -Bound Complex: A Small-Angle X-ray and Neutron Scattering Study." *Biochemistry* 33:8233-8239 (1994).
- Pearlstone and Smillie. "The interaction of rabbit skeletal muscle troponin-T fragments with troponin-I." *Canadian Journal of Biochemistry and Cell Biology* 63(3):212-218 (1985).
- Pearlstone and Smillie. "Effects of Troponin-I Plus -C on the Binding of Troponin-T and its Fragments to α -Tropomyosin." *J. Biol. Chem.* 258:2534-2542 (1983).
- Potter et al.. "Preparation of Troponin and Its Subunits." *Methods in Enzymology* 85:241-263 (1982).
- Pluskal et al.. "Immobilon™ PVDF Transfer Membrane: A New Membrane Substrate For Western Blotting of Proteins." *BioTechniques* 4(3):272-283 (1986).
- Schreier et al.. "Cloning, Structural Analysis, and Expression of the Human Slow Twitch Skeletal Muscle/Cardiac Troponin C Gene." *J. Biol. Chem.* 265:21247-21253 (1990).
- Stull and Buss. "Phosphorylation of Cardiac Troponin by Cyclic Adenosine 3':5'-Monophosphate-dependent Protein Kinase." *J. Biol. Chem.* 252:851-857 (1977).
- Townsend et al.. "Human Cardiac Troponin T: Identification of Fetal Isoforms and Assignment of the TNNT2 Locus to Chromosome 1q." *Genomics* 21:311-316 (1994).
- Vallins et al.. "Molecular cloning of human cardiac troponin I using polymerase chain reaction." *FEBS Letters* 270(1, 2):57-61 (1990).
- Wilkinson and Grand. "Comparison of amino acid sequence of troponin I from different striated muscles." *Nature* 271:31-35 (1978).
- Wu et al.. "Creatine Kinase MB Isoforms in Patients with Skeletal Muscle Injury: Ramifications for Early Detection of Acute Myocardial Infarction." *Clin. Chem.* 38(12):2396-2400 (1992).
- Zabel et al.. "Analysis of Creatine Kinase, CK-MB, Myoglobin, and Troponin T Time-Activity Curves for Early Assessment of Coronary Artery Reperfusion After Intravenous Thrombolysis." *Circulation* 887:1543-1550 (1993).
- Zot and Potter. "Structural Aspects of Troponin-Tropomyosin Regulation of Skeletal Muscle Contraction." *Ann. Rev. Biophys. Biophys. Chem.* 16:535-539 (1987).
- Cummins et al.. "Cardiac Specific Troponin-I Radioimmunoassay in the Diagnosis of Acute Myocardial Infarction." *American Heart Journal* 113:1333-1344 (1987).
- Bodor et. al. *Clin Chem* 38(11), 2203-2214, 1992.
- Katus et. al. *Clin Chem* 38(3), 386-393, 1992.
- Larve et. al. *Mol. Immun.* 29(2), 271-278, 1992.
- Larue et. al. *Clin Chem.* 39(6), 972-979, 1993.
- Takahashi et. al. *Clin Chem* 39(6), 1133, 1993.
- Sigma Catalog, 1992.
- Tietz. *Text Book of Clinical Chemistry*, 1986, pp. 1810, 1816.
- Syska. *Febs Letters*, 1974, pp. 253-257, vol. 40:2.

Figure 1a

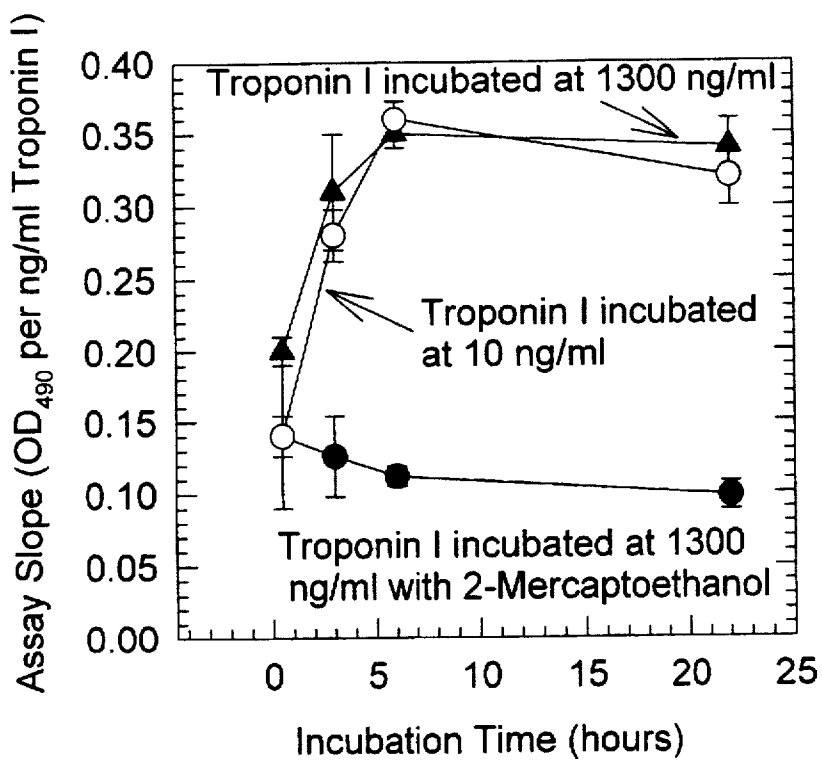


Figure 1b

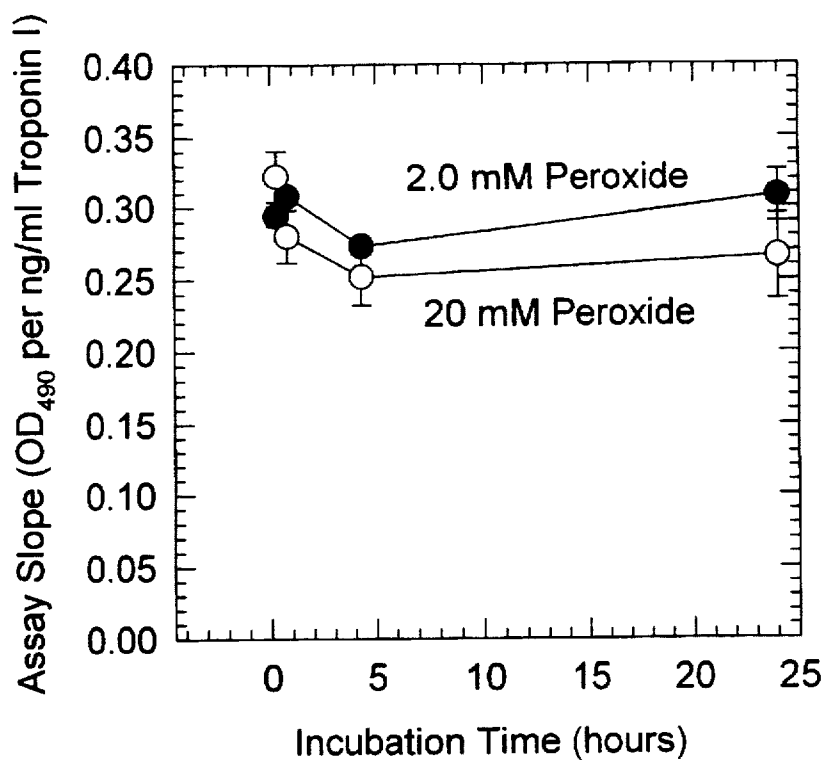


Figure 2

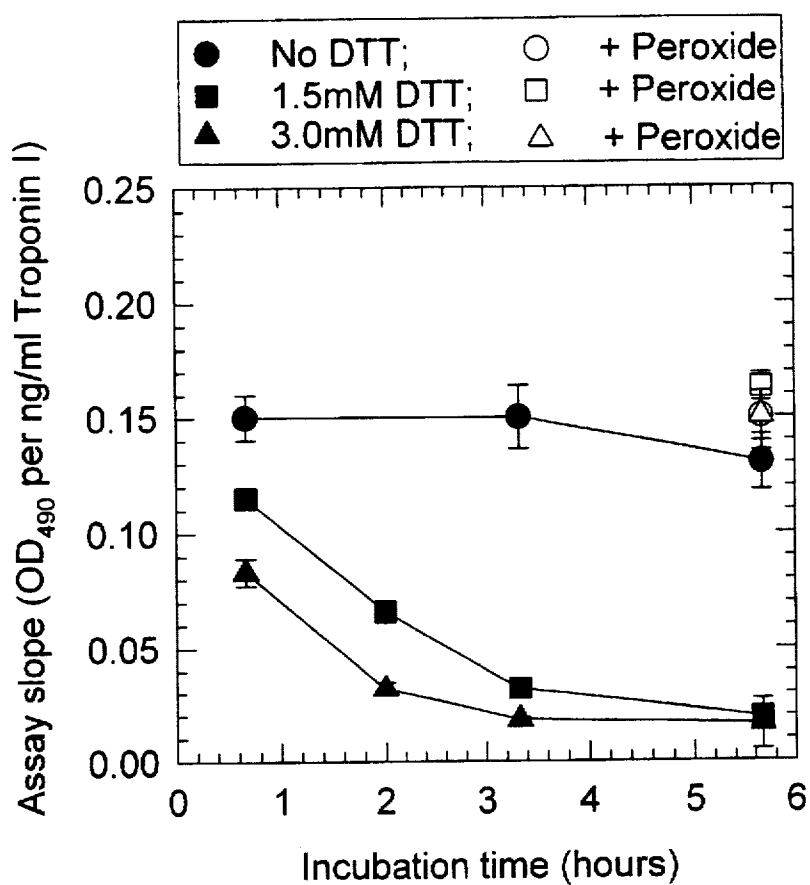


Figure 3

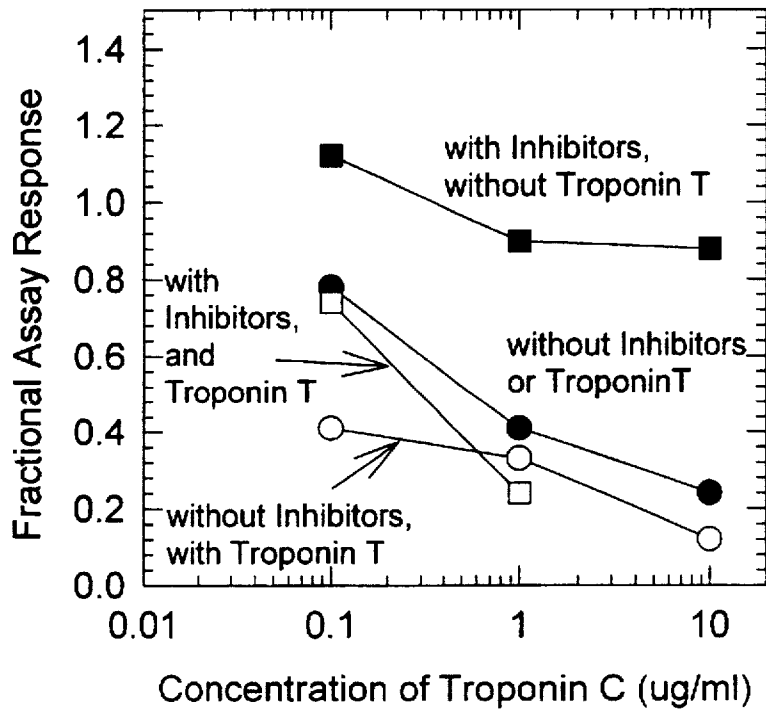
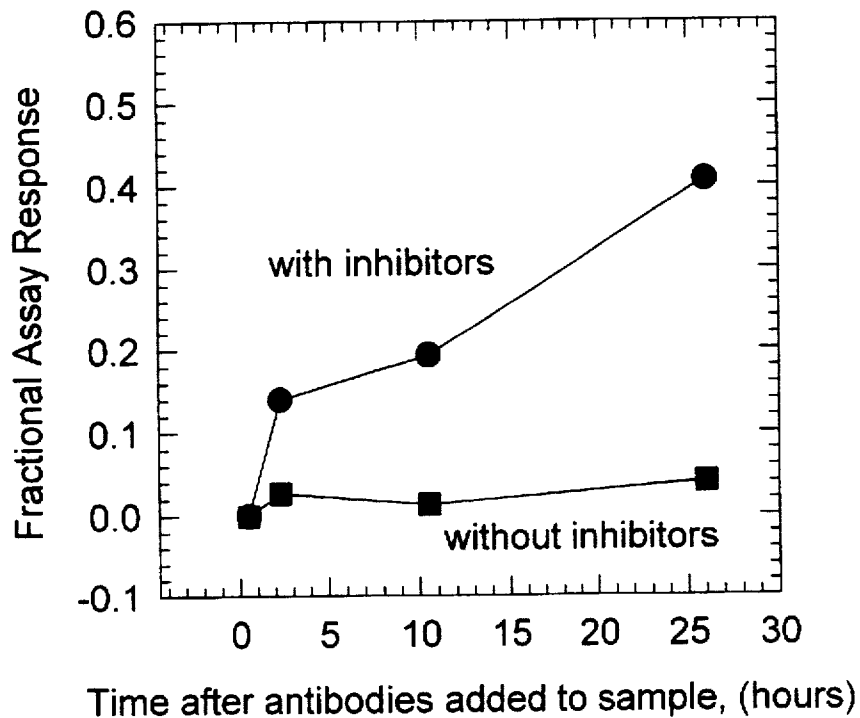
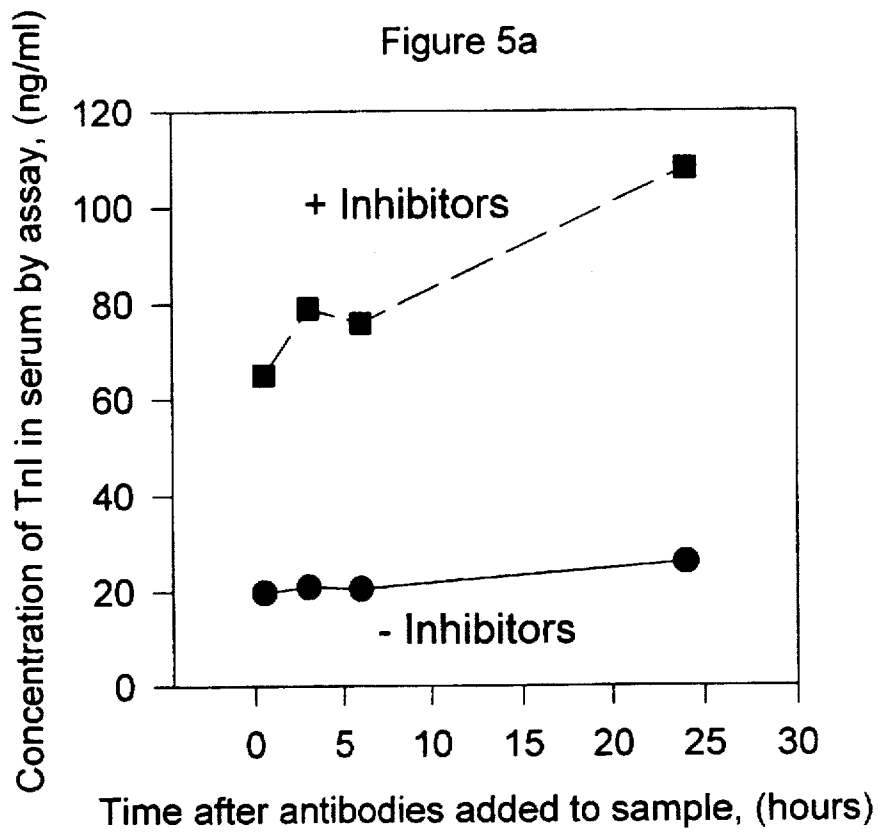
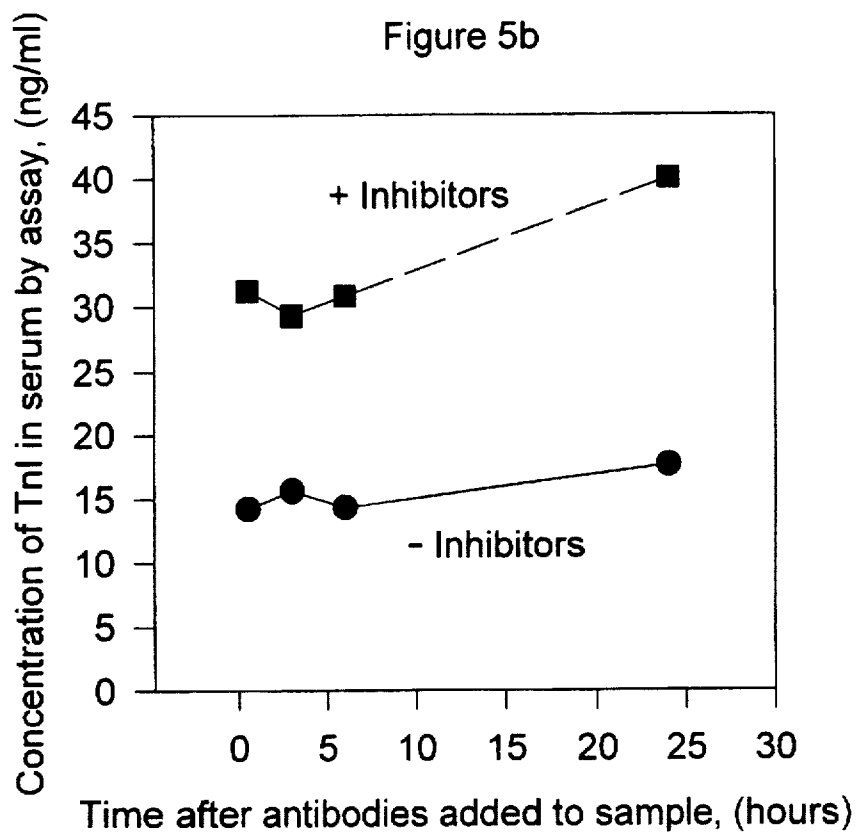


Figure 4







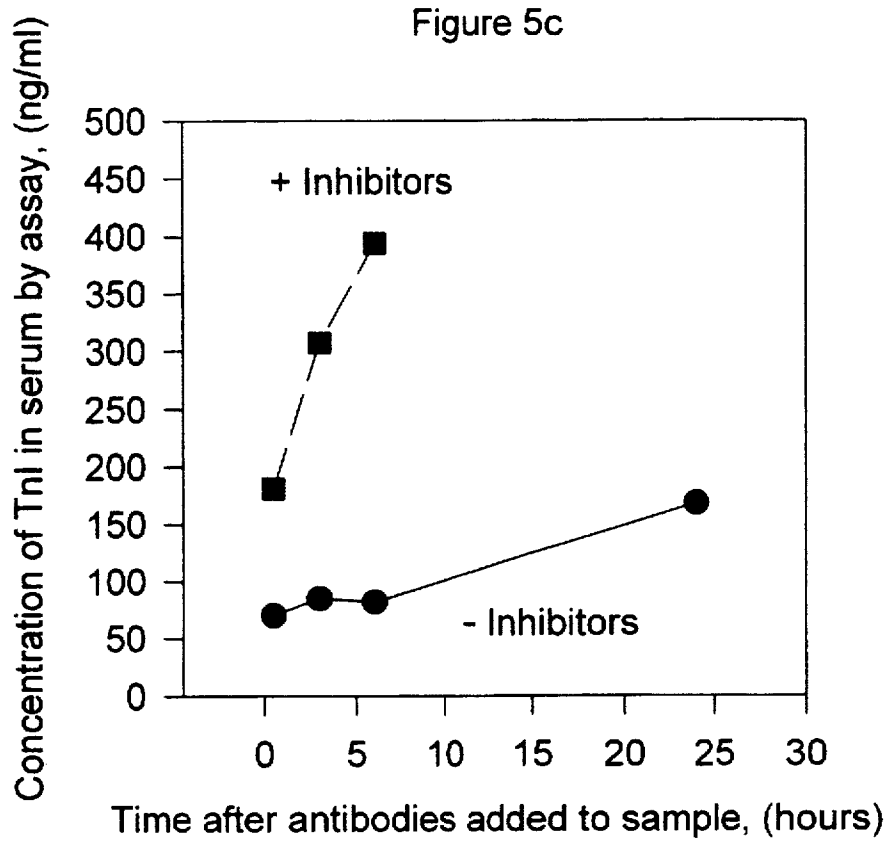


Figure 5d

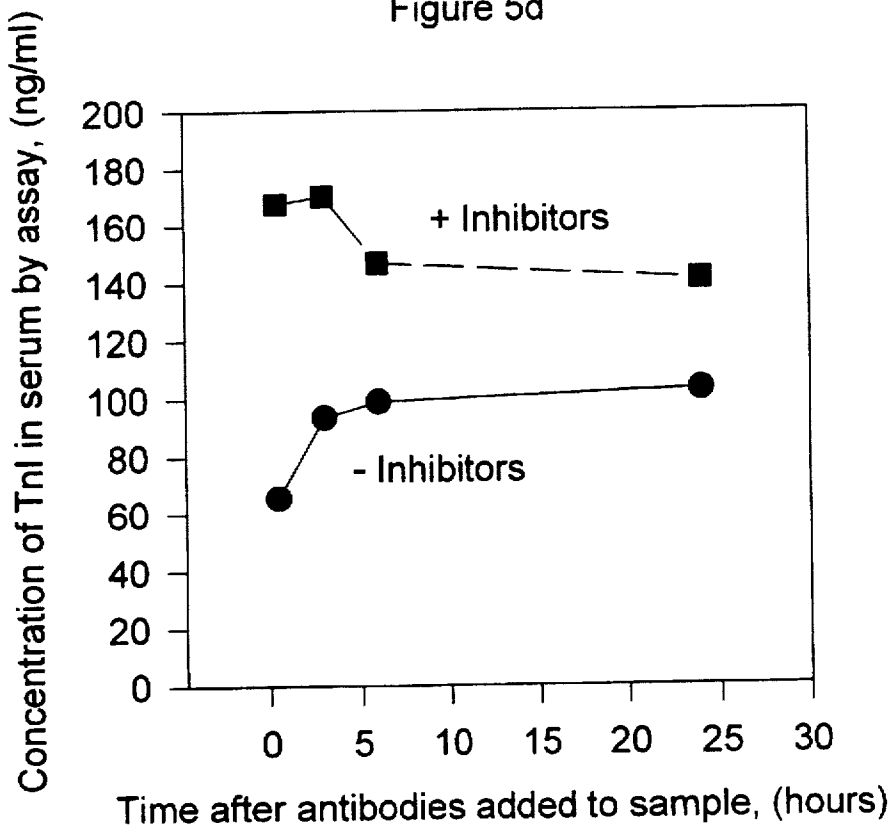


Figure 5e


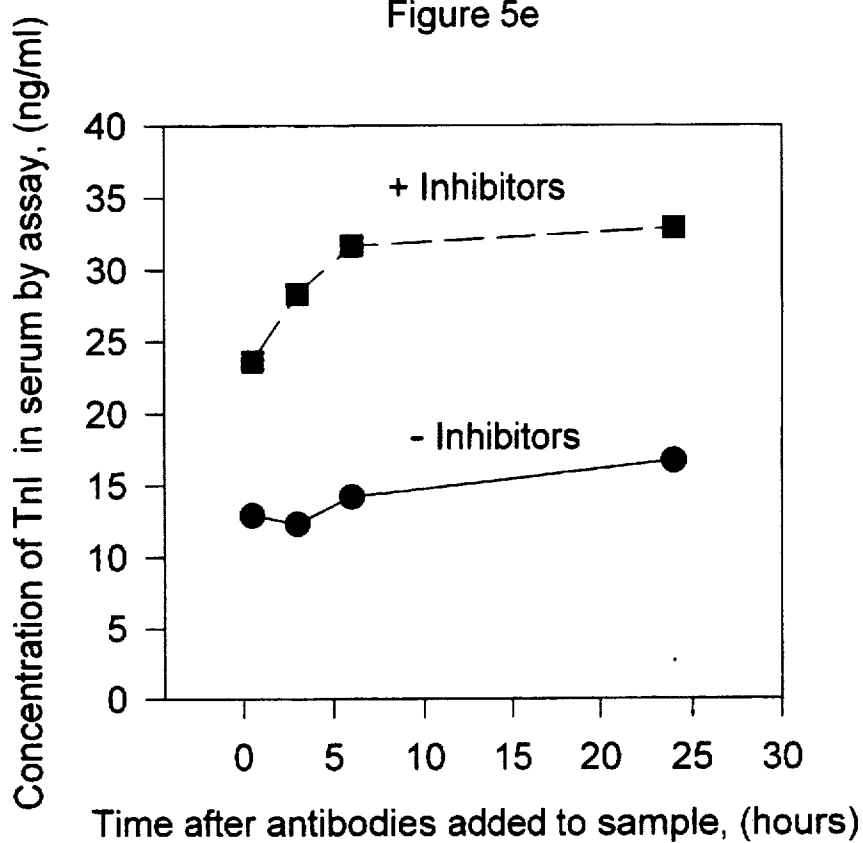


Figure 5f

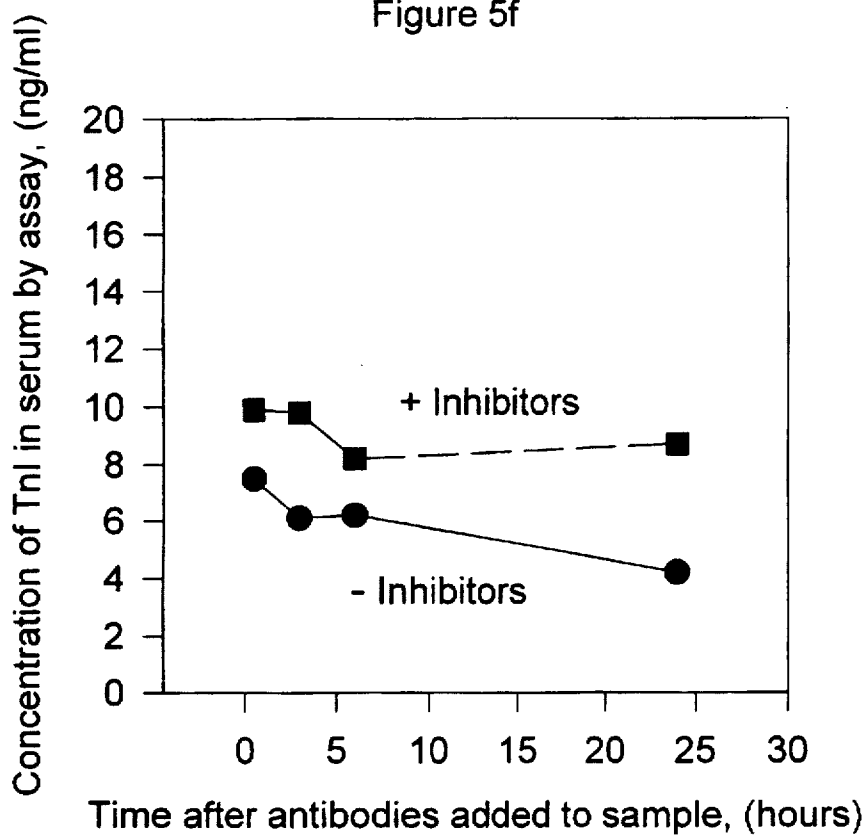

METHODS FOR THE ASSAY OF TROPONIN I AND T AND SELECTION OF ANTIBODIES FOR USE IN IMMUNOASSAYS

FIELD OF THE INVENTION

This invention relates to the assay of troponin I and troponin T in blood and more specifically to the changes in conformation of the proteins in blood and to the selection of antibodies to the various forms of the proteins and their use in immunoassays. In another aspect of the invention, compositions are taught for the stabilization and recovery of troponin I and T in immunoassays.

BACKGROUND

Myocardial infarction is one of the leading causes of death in the United States. Approximately 5 million individuals experiencing chest pain are evaluated every year in hospitals throughout the United States and of these people, less than 30% of them are subsequently found to have had a myocardial infarction. The accurate and rapid diagnosis of myocardial infarction is important both for the patient suffering a myocardial infarction and for the health care system which can minimize the costs incurred by rapidly identifying individuals who do need treatment.

The diagnosis of myocardial infarction is usually performed in the emergency department of a hospital. An individual having the symptoms of myocardial infarction is treated in different ways depending on the obviousness of the condition. Generally, an electrocardiogram is given to assess the condition of the heart; however, approximately 50% of patients experiencing myocardial infarction have a non-diagnostic electrocardiogram. The physician is then faced with a problem of diagnosing and treating the patient whom he expects may have a myocardial infarction.

The World Health Organization (WHO) has instituted guidelines for diagnosing myocardial infarction which state that an individual must exhibit 2 of the 3 following criteria: 1) chest pain or have a history of cardiac disease, 2) diagnostic electrocardiogram and 3) elevated creatine kinase (CK) or creatine kinase, MB isoenzyme (CKMB). For the 50% of the individuals who are presented to hospitals for a suspected myocardial infarction and who have a non-diagnostic electrocardiogram, the physician must rely on symptoms of chest pain and an elevated CK or CKMB to diagnose a myocardial infarction.

The assay of CK or CKMB is generally performed in hospital laboratories using sophisticated instrumentation. The assays include enzyme assays and immunoassays which detect the activity or mass of CK or CKMB present in blood samples.

During an episode of myocardial infarction, heart muscle cells die and release their contents to the blood stream. The CKMB, for example, which is released, among other cellular components, becomes elevated above a nominal value and can be diagnostic for myocardial infarction. The specificity of CKMB for diagnosing myocardial infarction is not 100% because another source of CKMB in the body is skeletal muscle. The mass of skeletal muscle in the body far exceeds the mass of cardiac muscle so that, through the normal catabolic turnover of skeletal muscle cells, the blood concentration of CKMB in healthy individuals will vary. In general, the concentration of CKMB which may be indicative of myocardial infarction is above 5-7 ng/ml (Circulation 87, 1542-1550 (1993), Clin. Chem. 39, 1725-1728 (1993)). The CKMB concentration of individuals who have skeletal muscle injury or who have exercised

has been reported to be elevated above 9 ng/ml (Clin. Chem. 38, 2396-2400 (1992)). The problem of specificity using CKMB as a marker for myocardial infarction has prompted the search for other specific markers which are released only from damaged heart muscle.

Troponin I and troponin T have recently been shown to be more specific than CKMB for diagnosing myocardial infarction (Circulation 83, 902-912 (1991), Clin. Chem. 40, 1291-1295 (1994)). Although troponin T has some disadvantages as a marker because it is elevated in patients experiencing renal disease (Clin. Chem. 41, 312-317 (1995)), the inventive methods herein successfully teach the use of troponin T as a diagnostic marker. The use of troponin I as a diagnostic marker for myocardial infarction also appears to meet many of the clinical requirements (Clin. Chem. 40, 1291-1295 (1994), Clin. Chem. 41, 312-317 (1995)).

The teachings of the instant invention provide methods for the selection of antibodies and their use in immunoassays for troponin I and troponin T. These proteins, along with troponin C, exist in both cardiac and skeletal muscle mainly as a ternary complex. In the muscle, the troponin complex is bound to tropomyosin which is, in turn, bound to the actin comprising the thin filaments. The state of troponin I and troponin T, whether free or bound as binary or ternary complexes, which are released from the muscle, has not been investigated.

The conformations of troponin I, T and C change upon binding when forming binary and ternary complexes (Biochemistry 33, 12800-12806 (1994), J. Biol. Chem. 254, 350-355 (1979), Ann. Rev. Biophys. Biophys. Chem. 16, 535-559 (1987)). An understanding of the conformational changes of troponin I and troponin T and the heterogeneity of the proteins in the blood is critical for the development of accurate diagnostic procedures for measuring troponin I and troponin T concentrations. In addition, troponin I is reported to be unstable in blood (Direction Insert for Troponin I Immunoassay, Sanofi/ERIA Diagnostics Pasteur, Marnes la Coquette, France), and the mechanisms responsible for the instability have not been understood. This invention addresses these problems and provides for stable troponin I and T compositions which are useful in immunoassays.

DESCRIPTION OF FIGURES

FIG. 1a illustrates the kinetics of air oxidation of troponin I as measured by immunoassay.

FIG. 1b illustrates the kinetics of oxidation by peroxide of troponin I as measured by immunoassay.

FIG. 2 illustrates the kinetics of reduction by dithiothreitol of troponin I and reoxidation of reduced troponin I by peroxide as measured by immunoassay.

FIG. 3 illustrates the effect of troponin C on the immunoassay of troponin I in the presence and absence of troponin T and binding inhibitors.

FIG. 4 illustrates the kinetics of disruption of human cardiac troponin ternary complex in the presence and absence of binding inhibitors as measured by an immunoassay for troponin I.

FIG. 5a illustrate the effect of binding inhibitors on troponin I immunoassays from a patient sample with confirmed myocardial infarction.

FIG. 5b illustrates the effect of binding inhibitors on troponin I immunoassays from a patient sample with confirmed myocardial infarction.

FIG. 5c illustrates the effect of binding inhibitors on troponin I immunoassays from a patient sample with confirmed myocardial infarction.

FIG. 5d illustrates the effect of binding inhibitors on troponin I immunoassays from a patient sample with confirmed myocardial infarction.

FIG. 5e illustrates the effect of binding inhibitors on troponin I immunoassays from a patient sample with confirmed myocardial infarction.

FIG. 5f illustrates the effect of binding inhibitors on troponin I immunoassays from a patient sample with confirmed myocardial infarction.

DETAILED DESCRIPTION

This invention is directed to the assay of troponin I and troponin T in body fluids, particularly, in human blood, serum and plasma. The presence of cardiac troponin I and T in the blood, above a nominal concentration, is diagnostic for damaged heart muscle. The teachings of this invention show that troponin I and T exist in various conformations in the blood which may be the same or different than their native conformations in muscle tissue. The ratios of the monomeric troponin I and T and the binary and ternary complexes may be related to the metabolic state of the heart. Based on the reactivities of antibodies to troponin I and T and to purified complexes of the troponins, the extent of troponin I and T complexation can now be elucidated in blood samples from patients suffering from myocardial infarction. The embodiments of this invention relate to the conformations of troponin I and T in blood and antibodies which recognize those conformations. Specifically, antibodies which recognize troponin I and T in the following forms are preferred: 1) The conformations of troponin I having intramolecularly oxidized and reduced cysteines, 2) The binary complexes of troponin I and T, of troponin I and C, of troponin T and C, and 3) The ternary complex of troponin I, T and C. In addition, methods are described for the improved recovery of troponin I and T in immunoassays. This invention answers an unmet need for the assays of troponin I and T in blood.

The troponin complex in muscle is comprised of troponin I, C and T and these components exist as various tissue specific isoforms. Troponin C exists as two isoforms, one from cardiac and slow-twitch muscle and one from fast-twitch muscle. Troponin I and T are expressed as different isoforms in slow-twitch, fast-twitch and cardiac muscle (Biochem. J. 171, 251-259 (1978), J. Biol. Chem. 265, 21247-21253 (1990), Hum. Genet. 88, 101-104 (1991), Circul. Res. 69, 1226-1233 (1991)). The unique cardiac isoforms of troponin I and T allow them to be distinguished immunologically from the other troponins of skeletal muscle. Therefore, the release into the blood of troponin I and T from damaged heart muscle has been related to cases of unstable angina and myocardial infarction. The prior art, however, has not addressed other forms of troponin I and T in blood.

The troponin complex in muscle is tightly bound to the contractile apparatus. Approximately 6% of the troponin T in cardiac tissue exists as an unbound protein in the cytoplasm and it is believed that this pool of troponin T is released from damaged muscle (Am. J. Cardiol. 67, 1360-1367 (1991)). The focus herein on troponin I and T for use as markers for myocardial infarction is based in part on size: because the proteins are relatively small, it is believed that they leak out of damaged cells faster than the larger proteins.

Antibodies to Troponin Complexes and to Troponin I and T

In a preferred embodiment, receptor proteins, for example, antibodies or binding fragments, are directed to the epitopes of troponin I which are insensitive to the oxidation state of the molecule. Human cardiac troponin I contains two cysteines, at positions 80 and 97 (FEBS Letters, 270, 57-61 (1990)). In the current art, during the purification of troponin I from tissues, the oxidation state of troponin I is directed toward the reduced form using various reductants, including mercaptoethanol, dithiothreitol and the like (Can. J. Biochem. 54, 546-553 (1976), Methods Enzymol. 85, 241-263 (1982)). After purification, the current art also teaches to maintain troponin I in the reduced form to prevent intermolecular disulfide formation (J. Biol. Chem. 258, 2951-2954 (1983)).

In the development of immunoassays for a target protein, the purified target protein acts as a standard with which to judge the sensitivity and specificity of the immunoassay using the antibodies which have been selected. The teachings of the instant invention show that the cysteines in troponin I can rapidly oxidize, intramolecularly, to alter the conformation of the protein. The degree of oxidation of troponin I has not previously been addressed with respect to its effect on the immunoassay process. The teachings described herein show that an apparent instability in the troponin I molecule is related to the dynamics of the intramolecular oxidation or reduction of the troponin I molecule. In addition, the selection of antibodies is taught for the accurate quantitation of troponin I in blood.

The terms sensitive and insensitive herein refer generally to the ability of an antibody to recognize or not to recognize a particular form of troponin I and T or complexes of troponin I and T. More specifically, a sensitive antibody which is useful in an immunoassay distinguishes one form or forms of troponin from another form and an insensitive antibody which is useful in an immunoassay does not distinguish one form or forms of troponin from another.

Purified, reduced troponin I undergoes an intramolecular oxidation of the cysteines, the rate of which is not dependent on the troponin I concentration. Special care should be exercised when preparing the oxidized troponin I form, especially in the presence of various thiol reducing agents, because of the possibility of forming mixed disulfides of the protein and the reducing thiol reagent. The mixed disulfide form of the protein may not behave as either the oxidized or reduced form of the molecule, especially if the antibodies used in the immunoassay bind to the region of the protein surrounding cysteines 80 and 97. Using purified preparations of oxidized and reduced troponin I, differential effects in the immunoassays using various antibodies raised to troponin I were observed. With some antibody pairs, the oxidized troponin I was hardly detectable in immunoassays, whereas with other pairs, the oxidation state had no effect on the immunoassay process. These results show that selection of antibodies to the troponin I molecule, without prior knowledge of the oxidation state of the troponin I, can result in antibodies and an immunoassay process which gives erroneous results. This conclusion is exemplified by immunoassays of troponin I from patients suffering myocardial infarction. The degree of oxidation of troponin I in patient samples is variable and indicates the possibility for apparent instabilities of troponin I assays of the current art.

In another preferred embodiment, antibodies or binding fragments are directed to the epitopes of the troponin I which are insensitive to the binding to troponin C and troponin T. The teachings herein show that antibodies are sensitive to the binding of troponin I to troponin C and to troponin T and

provide methods for the estimation of the extent of binding. Using purified preparations of troponin I, T and C, the effects of troponin I binding to troponin C and T on the recognition by antibody pairs is taught and is related to the dynamic state of troponin I in blood.

In another preferred embodiment, antibodies or binding fragments are directed to the epitopes of the troponin T which are insensitive of the binding to troponin C and troponin I. The teachings herein show that antibodies are sensitive to the binding of troponin T to troponin C and to troponin I and provide methods for the estimation of the extent of binding. Using purified preparations of troponin I, T and C, the effect of troponin T binding to troponin C and I on the recognition by antibody pairs is taught and is related to the dynamic state of troponin T in blood.

The degree of binding of troponin I and T to the components of the troponin complex or to other proteins of the contractile apparatus, including tropomyosin and actin, in blood can also be problematic for immunoassays depending on the degree and affinity of binding. In their native forms, the troponin complex exists in cardiac muscle and in slow- and fast-twitch skeletal muscle as a ternary complex of troponin I, C and T. Troponin I and T from skeletal muscle have different amino acid sequences than troponin I and T from cardiac muscle, respectively; however, troponin C from slow-twitch muscle has the same amino acid sequence as the cardiac muscle protein (Nature 271, 31-35 (1978), Arch. Biochem. Biophys. 186, 411-415 (1978), FEBS Lett. 292, 5-8 (1991)). The fast-twitch skeletal muscle troponin C, although not identical to the cardiac troponin C, can bind to cardiac troponin I (Biochemistry 33, 8464-8471 (1994), Proc. Natl. Acad. Sci. USA 90, 9036-9040 (1993)).

The release of troponin components, that is, troponin I, C and T, or components from the contractile apparatus, for example, tropomyosin and actin, from skeletal muscle, due to the normal turnover of skeletal muscle cells, may result in a significant amount of troponin and contractile apparatus components in the blood. Since skeletal muscle mass is much greater than cardiac muscle mass, the troponin components present in the blood of a normal individual may be derived largely from skeletal muscle. The circulating troponin components which are mainly derived from skeletal muscle would bind to cardiac troponin I and T which are released into the blood during a myocardial infarction or events which lead up to creating damaged heart muscle. As muscle damage progresses in an individual the troponin components derived from heart tissue will presumably rise in the blood. Thus, the concentration of troponin components in the blood from individuals experiencing a myocardial infarction may be differentially derived from both cardiac and skeletal muscle.

The form of troponin released from the heart, whether uncomplexed or as binary or ternary complexes, into the blood may indicate a particular condition of the heart. The assays taught herein provide for the analysis of release patterns which may allow the physician to diagnose a specific heart failure, for example, unstable angina as compared to myocardial infarction.

The clinical impact of an immunoassay measuring only the free troponin I or T from a patient experiencing a myocardial infarction can be very significant. Since the binding of troponin I and T to troponin components in the blood will be variable, depending on the troponin component concentrations, an analysis of the bound and free form of the troponin I and T in the blood must be considered. For example, the binding affinity of troponin I to troponin C, in the presence of calcium (which also is present in blood) is

$1.27 \times 10^8 \text{ M}^{-1}$ (Biochemistry 33, 12729-12734 (1994)). This implies that if the troponin C concentration is 100 ng/ml and the total (bound and free) troponin I concentration is 8 ng/ml, then the free concentration of troponin I is calculated to be 4.6 ng/ml. If the concentration of troponin I which is indicative of a myocardial infarction is 5 ng/ml or greater, then an assay which measures only the free form of troponin I, in this case, 4.6 ng/ml, will indicate to the physician that a myocardial infarction has not taken place. Generally, in hospital emergency departments which admit patients believed to have had a myocardial infarction, a blood sample from the individual will be obtained again in an hour or two if the first result is negative. In this example, the patient, having a total troponin I concentration of 8 ng/ml, (which is defined as positive for a myocardial infarction), but only a measured concentration of 4.6 ng/ml, (would be defined as a negative result), would not be treated and would continue to accrue damaged heart muscle. Interpretation of results of troponin T assays would also suffer from troponin T binding to components of the contractile apparatus in blood. Thus, immunoassays of the current art which measure the free troponin I and T may not correctly diagnose a myocardial infarction when the troponin I or T concentration, respectively, is near the decision point.

In a particularly preferred embodiment, antibodies or binding fragments are directed to the cardiac troponin complex. Specifically, antibodies are directed to cardiac specific epitopes of troponin I and T of the troponin complex or of the troponin I/T, I/C and T/C interfaces in the complex. The teachings herein show that antibodies which are raised to troponin I and T bind poorly to troponin I or T of the ternary complex. Furthermore, the teachings herein show that the troponin complex exists in the blood of patients who have experienced myocardial infarction. Methods are also described which teach one skilled in the art to assess the amount of troponin complex in the blood relative to the free troponin I and T or binary complexes of troponin I and T, using antibodies which bind to the free troponin molecules.

The association constants for the binary and ternary complexes of troponin will be altered during the immunoassay process because of the binding of antibodies to the troponin components and complexes. The change in the affinity constants during the immunoassay may be significant or insignificant and will be a function of the antibody concentrations, the affinity of the antibodies for the troponin components and complexes, and the time that the antibodies are allowed to bind to the troponin. These variables can change the perceived concentration of troponin I and T and lead to erroneous conclusions about the troponin concentration. For example, if two immunoassays utilize different antibody pairs for performing a sandwich immunoassay and their antibody concentrations and affinities for troponin I or T are different, and if a proportion of the troponin I or T occurs in the sample as binary and ternary complexes, one may expect that each immunoassay will give a different result. In addition, if blood samples contain varying concentrations of troponin C, then the proportion of troponin I and T that is bound to troponin C as a binary complex will differentially perturb each immunoassay.

The teachings of the instant invention demonstrate that the troponin ternary complex is more stable to dissociation than the binary complexes of troponin.

In another preferred embodiment, antibodies or binding fragments are directed to epitopes which are insensitive to proteolytic degradation of the N-terminal region of troponin I.

The conformation of troponin I is also reported to be affected by phosphorylation/dephosphorylation (Biophys. J.

63, 986-995 (1992), *Biochem.* 33, 12729-12734 (1994)). In another preferred embodiment, antibodies or binding fragments are directed to epitopes of either troponin I or troponin I complexes, which are sensitive or insensitive to the phosphorylation state of troponin I. Troponin I can be phosphorylated using methods described in *J. Biol. Chem.* 252, 851-857 (1977). The phosphorylated and dephosphorylated preparations of troponin I can be utilized as immunogens for generating antibodies as well as antigens for the selection of antibodies to the phosphorylated and dephosphorylated troponin I.

The troponin complex can be dissociated into the component proteins using various treatments, including high concentrations of urea, low pH and metal chelating agents which bind divalent metal cations, particularly calcium and magnesium (*Methods Enzymol.* 85, 241-263 (1982)). These treatments are, in general, very harsh and require several hours. Thus, these conditions for dissociating the troponin complex are not practical for immunoassays which must be performed in a matter of minutes on samples from individuals who may be suffering a myocardial infarction.

The generation and selection of antibodies which are preferentially either sensitive or insensitive to the binding of troponin I or T in binary complexes are accomplished by first preparing binary troponin I/T, T/C and I/C complexes from purified components (*J. Biol. Chem.* 254, 350-355 (1979), *J. Biol. Chem.* 258, 2534-2542 (1983), *J. Biol. Chem.* 258, 2951-2954 (1983), *Can. J. Biochem. Cell Biol.* 63, 212-218 (1985), *Biochemistry* 33, 12729-12734 (1994), *Ann. Rev. Biophys. Biophys. Chem.* 16, 535-559 (1987)). The generation and selection of antibodies which are sensitive or insensitive to the binding of troponin I or T in the ternary complex are accomplished, also, by first purifying the ternary complex (*Methods Enzymol.* 85, 241-263 (1983)) or by reconstitution of the complex using the purified troponin components. One skilled in the art will recognize that various other contractile apparatus proteins which may be associated with the binary or ternary complexes of troponins can also be constructed from the purified components and that the resultant complex can be utilized to generate and select antibodies as taught by the instant invention. The purified complexes are then injected, for example, into mice or rabbits, to generate monoclonal or polyclonal antibodies. One skilled in the art will recognize that many procedures are available for the production of antibodies, for example, as described in *Antibodies, A Laboratory Manual*, Ed Harlow and David Lane, Cold Spring Harbor Laboratory (1988), Cold Spring Harbor, N.Y. One skilled in the art will also appreciate that binding fragments which mimic antibodies can also be prepared from genetic information by various procedures (*Antibody Engineering: A Practical Approach* (Borrebaeck, C., ed.), 1995, Oxford University Press, Oxford; *J. Immunol.* 149, 9143-3920 (1992)).

The antibodies which are generated are selected by first screening for affinity and specificity with the purified binary or ternary complexes and comparing the results to the affinity and specificity of the antibodies with the purified troponin I and T molecules for the desired properties which are defined by the immunoassay process.

The screening procedure can involve immobilization of the purified troponin I or T or binary or ternary complexes or peptides corresponding to cardiac specific sequences of the troponins in separate wells of microtiter plates. The solution containing a potential antibody or groups of antibodies is then placed into the respective microtiter wells and incubated for about 30 min to 2 h. If an antibody to the

protein of interest is present in the solution, it will bind to the immobilized troponin. In screening antibodies for binding to interfaces of binary or ternary complexes of troponin, an antibody is first selected which binds to the binary or ternary complex immobilized in the microtiter well. That antibody is then further screened for its ability to bind to free troponin components; that is, the potential interface antibody should not bind to free troponin I, C or T which is immobilized in microtiter wells. In addition, the interface antibody should never be capable of forming a sandwich assay with binary or ternary complexes and an antibody which is known to bind to a specific troponin component in the complex in the presence of binding inhibitors which are known to disrupt the troponin complex. If this latter condition is met, then the potential interface antibody should also not be capable of forming a sandwich assay with the troponin complex and an antibody to a different troponin component than was used in the previous screen in the presence of binding inhibitors. If this condition is also met, then an interface antibody has been selected for a binary complex. An extra immunoassay must be performed for selecting an interface antibody to a ternary complex; that is, if the previous two conditions are met, then the potential interface antibody should also not be capable of forming a sandwich assay with the troponin complex and an antibody to a different troponin component than was used in the two previous screens. The microtiter wells are then washed and a labeled secondary antibody (for example, an anti-mouse antibody conjugated to alkaline phosphatase if the raised antibodies are mouse antibodies) is added to the wells and incubated for about 30 min and then washed. Substrate is added to the wells and a color reaction will appear where antibody to the troponin is present. The antibodies which are of interest are then further analyzed for affinity and specificity to the cardiac specific molecules and for complementarity in forming sandwich complexes with the antigens. Those skilled in the art will recognize that many approaches can be taken in producing antibodies or binding fragments and screening and selecting for affinity and specificity for the various troponin antigens, but these approaches do not change the scope of the invention.

Assays for Troponin Complexes and Uncomplexed Troponin I and T

A particularly preferred embodiment of this invention is directed to the assay of troponin I and troponin T, particularly immunoassays, wherein the antibodies selected for the assay bind to cardiac specific sequences of the ternary complex, of the binary complexes and of the uncomplexed (free) troponin I or T in order to measure the complexed (bound) and free fractions of troponin I and T, respectively. The cardiac specific sequences of troponin I and T are described in *FEBS Lett.* 270, 57-61 (1990) and *Genomics* 21, 311-316 (1994). A synthetic peptide comprised of 14 amino acids which mimics a cardiac specific sequence of troponin I and methods used to prepare antibodies to the peptide are described in an International Patent Application number PCT/US94/05468 International Publication Number WO 94/27156.

The immunoassay can be formulated with a cocktail of antibodies to bind all the troponin complexes and the unbound troponin I and T. Alternatively, the immunoassay can be formulated with antibodies that recognize epitopes of the troponin I and T in the complexes and also the unbound troponin I and T. In addition, the immunoassay can be formulated with antibodies that bind epitopes at interfaces of the component proteins in the complexes and antibodies that bind the unbound troponin I and T.

A preferred immunoassay for troponin I or T involves conjugation of a cocktail of antibodies to a label or a signal

generator to form antibody conjugates, which are capable of binding to cardiac specific regions of the troponin complexes of troponin I or T and to unbound troponin I or T. Another antibody or cocktail of antibodies is immobilized on a solid phase, for example, a membrane as taught in BioTechniques 4, 272-283 (1986), and the membrane is placed in a device, for example, as described in allowed U.S. patent, Ser. No. 07/961,267, now pending as Ser. No. 08/458,276. The immobilized antibody is complementary to the antibody conjugate and forms sandwich complexes with the troponin I or T complexes and the troponin I or T, respectively. A plasma or serum sample suspected of containing troponin complexes or components from damaged heart muscle is mixed with the antibody conjugate to form a reaction mixture and it is allowed to incubate. The reaction mixture is then applied to the aforementioned device. The sample flows through the membrane and the troponin complexes and components, bound to the antibody conjugates, bind to the immobilized antibodies and excess, unbound antibody conjugate is washed away with a wash buffer. The signal is developed and read, either visually or instrumentally.

A particularly preferred immunoassay for troponin I involves conjugation of at least two antibodies to a label or a signal generator to form an antibody conjugate. One of the conjugate antibodies is capable of binding to the troponin T component of the troponin complexes and the other antibody is capable of binding to the free and binary troponin I molecules. Another antibody or cocktail of antibodies is immobilized on a solid phase, for example, a membrane, and the membrane is placed in a device, as described previously. The immobilized antibody is complementary with the antibody conjugate antibodies to form sandwich complexes with either troponin I bound to troponin complexes or to the uncomplexed troponin I. A plasma or serum sample suspected of containing troponin complexes or components from damaged heart muscle is mixed with the antibody conjugate to form a reaction mixture and it is allowed to incubate. The reaction mixture is then applied to the aforementioned device. The sample flows through the membrane and the troponin complexes and components, bound to the antibody conjugates, bind to the immobilized antibodies and excess, unbound antibody conjugate is washed away with a wash buffer. The signal is developed and read, either visually or instrumentally. In this assay procedure, the antibody conjugate binds to the troponin complexes through the troponin T specific antibody and all free and binary troponin I molecules through the troponin I specific antibody. The capture antibody or antibodies on the solid phase bind antibody conjugates having bound only free troponin I or troponin complexes containing troponin I.

Another particularly preferred immunoassay for troponin I involves conjugation of an antibody or a cocktail of antibodies to a label or a signal generator to form an antibody conjugate. The antibody conjugate binds to either troponin I bound to troponin complexes or to the uncomplexed troponin I. Immobilized on a solid phase, for example, a membrane, in 3 discrete zones, are antibodies or cocktails of antibodies which bind the ternary complex, the binary complexes of troponin I and the free troponin I, and the membrane is placed in a device, as described previously. For example, a troponin T antibody which binds to the troponin T of the ternary complex is immobilized in one discrete zone, a troponin I antibody which binds to the troponin I binary complexes (troponin I/C and I/T) is immobilized in another discrete zone and a troponin I antibody which binds to only the uncomplexed troponin I is

immobilized in yet another discrete zone. The immobilized antibodies are complementary with the antibody conjugate antibodies to form sandwich complexes with complexed or uncomplexed troponin I, as defined by each discrete zone. Alternatively, immobilized on a solid phase, for example, a membrane, in 2 discrete zones, are antibodies or cocktails of antibodies which bind the troponin I complexes (binary and ternary) and the free troponin I. For example, a troponin T antibody which binds to the troponin T of the ternary complex and a troponin I antibody which binds to the troponin I of the binary complexes are immobilized in one discrete zone and a troponin I antibody which binds to only the uncomplexed troponin I is immobilized in the second discrete zone. The immobilized antibodies are complementary with the antibody conjugate antibodies to form sandwich complexes with complexed or uncomplexed troponin I, as defined by each discrete zone. A further embodiment of this invention utilizes antibodies on the solid phase for detection of troponin I complexes which bind to the interfaces of the binding domains of troponin I/T and I/C. A plasma or serum sample suspected of containing troponin complexes or components from damaged heart muscle is mixed with the antibody conjugate to form a reaction mixture and it is allowed to incubate. The reaction mixture is then applied to the aforementioned device. The sample flows through the membrane and the troponin complexes and components, bound to the antibody conjugate(s), bind to the respective immobilized antibodies in the discrete zones and excess, unbound antibody conjugate is washed away with a wash buffer. The signal is developed and read, either visually or instrumentally. In this assay procedure, the antibody conjugate binds to the troponin I and the troponin I binary and ternary complexes through the troponin I specific antibody or antibodies. The capture antibody or antibodies in discrete zones on the solid phase bind the antibody conjugates which are specific to the uncomplexed troponin I or troponin complexes containing troponin I. This immunoassay allows quantification of the fractions of troponin I, namely, the complexed and the uncomplexed fractions. The inventive teachings described herein show that uncomplexed and complexed troponin exists in plasma and serum samples from patients with confirmed myocardial infarction. The determination of the complexed and uncomplexed troponin I fractions may yield important clinical data relating to the type and extent of muscle damage, for example, from unstable angina or myocardial infarction or to the success of thrombolytic therapy.

A particularly preferred immunoassay for troponin T involves conjugation of at least two antibodies to a label or a signal generator to form an antibody conjugate. One of the conjugate antibodies is capable of binding to the troponin I component of the troponin complexes and the other antibody is capable of binding to the free and binary troponin T molecules. Another antibody or cocktail of antibodies is immobilized on a solid phase, for example, a membrane, and the membrane is placed in a device, as described previously. The immobilized antibody is complementary with the antibody conjugate antibodies to form sandwich complexes with either troponin T bound to troponin complexes or to the uncomplexed troponin T. A plasma or serum sample suspected of containing troponin complexes or components from damaged heart muscle is mixed with the antibody conjugate to form a reaction mixture and it is allowed to incubate. The reaction mixture is then applied to the aforementioned device. The sample flows through the membrane and the troponin complexes and components, bound to the antibody conjugates, bind to the immobilized antibodies and

excess, unbound antibody conjugate is washed away with a wash buffer. The signal is developed and read, either visually or instrumentally. In this assay procedure, the antibody conjugate binds to the troponin complexes through the troponin I specific antibody and all free and binary troponin T molecules through the troponin T specific antibody. The capture antibody or antibodies on the solid phase bind antibody conjugates having bound only free troponin T or troponin complexes containing troponin T.

Another particularly preferred immunoassay for troponin T involves conjugation of an antibody or a cocktail of antibodies to a label or a signal generator to form an antibody conjugate. The antibody conjugate binds to either troponin T bound to troponin complexes or to the uncomplexed troponin T. Immobilized on a solid phase, for example, a membrane, in 3 discrete zones, are antibodies or cocktails of antibodies which bind the ternary complex, the binary complexes of troponin T and the free troponin T, and the membrane is placed in a device, as described previously. For example, a troponin I antibody which binds to the troponin I of the ternary complex is immobilized in one discrete zone, a troponin T antibody which binds to the troponin T binary complexes (troponin I/T and C/T) is immobilized in another discrete zone and a troponin T antibody which binds to only the uncomplexed troponin T is immobilized in yet another discrete zone. The immobilized antibodies are complementary with the antibody conjugate antibodies to form sandwich complexes with complexed or uncomplexed troponin T, as defined by each discrete zone. Alternatively, immobilized on a solid phase, for example, a membrane, in 2 discrete zones, are antibodies or cocktails of antibodies which bind the troponin T complexes (binary and ternary) and the free troponin T. For example, a troponin I antibody which binds to the troponin I of the ternary complex and a troponin T antibody which binds to the troponin T of the binary complexes are immobilized in one discrete zone and a troponin T antibody which binds to only the uncomplexed troponin T is immobilized in the second discrete zone. The immobilized antibodies are complementary with the antibody conjugate antibodies to form sandwich complexes with complexed or uncomplexed troponin T, as defined by each discrete zone. A further embodiment of this invention utilizes antibodies on the solid phase for detection of troponin T complexes which bind to the interfaces of the binding domains of troponin C/T and I/T. A plasma or serum sample suspected of containing troponin complexes or components from damaged heart muscle is mixed with the antibody conjugate to form a reaction mixture and it is allowed to incubate. The reaction mixture is then applied to the aforementioned device. The sample flows through the membrane and the troponin complexes and components, bound to the antibody conjugate(s), bind to the respective immobilized antibodies in the discrete zones and excess, unbound antibody conjugate is washed away with a wash buffer. The signal is developed and read, either visually or instrumentally. In this assay procedure, the antibody conjugate binds to the troponin T and the troponin T binary and ternary complexes through the troponin T specific antibody or antibodies. The capture antibody or antibodies in discrete zones on the solid phase bind the antibody conjugates which are specific to the uncomplexed troponin T or troponin complexes containing troponin T. This immunoassay allows quantification of the fractions of troponin T, namely, the complexed and the uncomplexed fractions. The inventive teachings described herein show that uncomplexed and complexed troponin exists in plasma and serum samples from patients with confirmed myocardial

infarction. The determination of the complexed and uncomplexed troponin T fractions may yield important clinical data relating to the type and extent of muscle damage, for example, from unstable angina or myocardial infarction or to the success of thrombolytic therapy.

Another particularly preferred immunoassay independently measures the cardiac troponin ternary complex, the cardiac troponin binary complexes (troponin I/T, T/C and I/C) and the free cardiac troponin I and T. This method involves conjugation of antibodies or a cocktail of antibodies to a label or a signal generator to form an antibody conjugate. The antibody conjugates bind to either troponin T and I bound to troponin complexes or to the uncomplexed troponin T and I. Immobilized on a solid phase, for example, a membrane, in 6 discrete zones, are antibodies or cocktails of antibodies which bind the ternary complex, the binary complexes of troponin I and T and the free troponin I and T, and the membrane is placed in a device, as described previously. For example, a troponin I antibody which binds to the troponin I of the ternary complex is immobilized in one discrete zone, 3 different antibodies, each recognizing the interfaces of the binding domains of troponin I/T, T/C and I/C, are immobilized in 3 discrete zones, a troponin I antibody which binds to the free troponin I is immobilized in another zone and a troponin T antibody which binds to the free troponin T is immobilized in another zone. The immobilized antibodies are complementary with the antibody conjugate antibodies to form sandwich complexes with complexed and uncomplexed troponin T and I, as defined by each discrete zone. A plasma or serum sample suspected of containing troponin complexes or components from damaged heart muscle is mixed with the antibody conjugate to form a reaction mixture and it is allowed to incubate. The reaction mixture is then applied to the aforementioned device. The sample flows through the membrane and the troponin complexes and components, bound to the antibody conjugate(s), bind to the respective immobilized antibodies in the discrete zones and excess, unbound antibody conjugate is washed away with a wash buffer. The signal is developed and read, either visually or instrumentally. In this assay procedure, the antibody conjugates bind to the free troponin I and T, to the troponin I and T of the binary complexes and the troponin I or T of the ternary complex. The capture antibodies in discrete zones on the solid phase bind the antibody conjugates which are specific to the free troponin I and T and to the troponin complexes. This immunoassay allows quantification of each of the ternary troponin complex, the troponin T/C, I/T, I/C and the free troponin I and T. The inventive teachings described herein show that uncomplexed and complexed troponin exists in plasma and serum samples from patients with confirmed myocardial infarction. The determinations of the individual troponin ternary complex, the troponin I and T binary complexes and uncomplexed troponin I and T fractions may yield important clinical data relating to the type and extent of muscle damage, for example, from unstable angina or myocardial infarction or to the success of thrombolytic therapy.

In another embodiment of this invention, inhibitors which affect the affinity constants of the association of troponin I complexes or of troponin T complexes are added to the sample prior to formation of the reaction mixture so that the free troponin I or troponin T is measured, respectively. The binding inhibitors may disrupt the troponin complexes or they may open up or partially unravel the complex, such that some or all interactions between the troponin components are broken so that epitopes may be more easily assessable to

the antibodies for binding. Inhibitors can be selected from the group of compounds which comprise, but is not limited to, metal chelating agents and peptides which compete with troponin I or troponin T for binding to proteins of the contractile apparatus. Metal chelating agents, particularly those which bind to calcium and magnesium, lower the affinity constant, for example, of troponin I and troponin C binding by about a factor of 10 as compared to the affinity in the presence of calcium (Biochemistry 33, 12729-12734 (1994)). Peptides which affect troponin I and troponin T binding to proteins of the contractile apparatus include mastoparan, melittin and peptide sequences which mimic the troponin I and T sequences at their binding domains with the proteins of the contractile apparatus ((Biochemistry 31, 11326-11334 (1992), J. Biol. Chem. 267, 15715-15723 (1992), Biochemistry 33, 8233-8239 (1994)). Other peptides which are useful as inhibitors are those which mimic the binding domains of the troponin components. The binding domains are described, for example, in Ann. Rev. Biophys. Biophys. Chem. 16, 535-559 (1987), and with the binding domain information, one skilled in the art synthesizes the peptide which mimics the peptide of the protein at the binding domain.

In another embodiment of this invention, troponin C is added to the sample prior to the formation of the reaction mixture of the immunoassay so that all or nearly all of the troponin I or T in the sample will be bound by troponin C during the course of the assay. The troponin C concentration in the sample should be about 0.5 $\mu\text{g/ml}$ to 100 $\mu\text{g/ml}$ and preferably about 1 to 10 $\mu\text{g/ml}$.

In another embodiment of this invention, troponin C and T are added to samples prior to formation of the reaction mixture of immunoassays for troponin I in order to bind all or nearly all of the troponin I in the form of the ternary troponin complex. The troponin C and T concentrations in the sample should be about 0.5 to 100 $\mu\text{g/ml}$ and preferably about 1 to 10 $\mu\text{g/ml}$.

In yet another embodiment of this invention, troponin C and I are added to samples prior to formation of the reaction mixture of immunoassays for troponin T in order bind all or nearly all of the troponin T in the form of the ternary troponin complex. The troponin C and T concentrations in the sample should be about 0.5 to 100 $\mu\text{g/ml}$ and preferably about 1 to 10 $\mu\text{g/ml}$.

These embodiments wherein troponin components are added to the sample prior to formation of the reaction mixture have several advantages.

Firstly, troponin I adsorbs tenaciously to glass surfaces and various membranes which can result in a lower measured troponin I concentration. Troponin T also adsorbs to surfaces. However, when bound to troponin C or in the ternary complex, the adsorptive characteristics of troponin I and T are dramatically reduced. Thus, the recovery of the troponin I/C or T/C complex or the ternary complex can be better than troponin I or T. In this embodiment, antibodies that recognize the troponin I or T complexes are used in the immunoassays.

Secondly, if antibodies which bind only to the complexed troponin I or T are required, then the antibody selection process is less stringent because the antibodies are not required to have a similar affinity to the free troponin I or T.

One skilled in the art will appreciate the inventive teachings described herein and will recognize with these teachings that addition of reagents to a device or to each other, as recited in the embodiments, has many forms, and these forms are within the scope of this invention.

Stabilization of Troponin or T for Calibrator Reagents

Troponin I and T are known to be unstable in aqueous formulations, as well as in patient samples. The (apparent) instabilities of the proteins, as taught herein, are related to the oxidation state of the troponin I, the propensity of troponin I and T to form complexes with other troponin proteins and the adsorptive characteristics of troponin I and T.

Stabilization of troponin I is performed by the intramolecular oxidation of the cysteines and the protein is stored without thiol reducing agents, such as mercaptoethanol, dithiothreitol and the like.

The storage of troponin I in solutions containing high concentrations of thiol reductants will maintain the cysteines, overall, in the reduced form. However, intramolecular oxidation and reduction of a protein is a dynamic process whereby the protein will exist for some time in the oxidized form even in the presence of the reductants. In the case of reductants, such as mercaptoethanol or N-acetylcysteine, that is, reductant molecules with only a single thiol group, mixed disulfides of the reductant and the protein thiol will form. The half-life of this mixed disulfide will be a function of the reductant concentration and the rate of intramolecular oxidation; that is, the mixed disulfide can be reduced by both the thiol reductant reagent and the other protein cysteine, assuming that both cysteines are not in the mixed disulfide form. In the case of reducing the intramolecularly oxidized troponin I, the reductant with a single thiol group will reduce the intramolecular cystine to yield a cysteine and a mixed disulfide of the protein. The mixed disulfide of the protein will be reduced by either the cysteine of the protein or the thiol reductant. This process continues and eventually results in depleting the reductant concentration to a level where it can no longer maintain the protein in the reduced state. As the reductant concentration approaches the concentration of thiol in the troponin I, the protein cysteine and the thiol reductant reagent can form mixed disulfides, which will not be reduced by the thiol reductant. Alternatively, the protein will oxidize, intramolecularly, and the thiol reductant is not in sufficient concentration to reduce the cystine. The end result, upon depletion of the thiol reductant, will be a mixture of troponin I which is in the intramolecularly oxidized form and protein which is in the mixed disulfide form. Each of these forms of troponin I has a different conformation.

In the case of utilizing thiol reductant reagents which possess two thiol groups, for example dithiothreitol or dithioerythritol, the end result, upon depletion of the thiol reductant, will be only the oxidized form of the troponin I.

Therefore, antibodies which are sensitive to the oxidation state of the troponin I will differentially recognize the various forms of the troponin I in the immunoassay. The immunoassay will then measure an inaccurate concentration of troponin I.

The teachings herein describing the oxidation of troponin I and the adsorptive behavior of troponin I and T have led to the finding that troponin I and T can be stabilized as the troponin ternary complex.

A preferred composition of stabilized troponin I comprises an aqueous solution of the intramolecularly oxidized troponin I.

A particularly preferred composition of stabilized troponin I and T comprises a buffered solution of the ternary complex of troponin I, T and C in the presence or absence of calcium and magnesium salts. A preferred range of pH of the solution is between 6 and 9 and a range of calcium and magnesium salts concentrations, for example calcium and

magnesium chloride, of between 0.01 mM and 10 mM. A particularly preferred buffered solution consists of up to about 100% human serum or plasma. The ternary complex can be formed from the component troponin I, T and C, or alternatively, it can be isolated and purified from cardiac or skeletal muscle (Methods Enzymol. 85, 241-263 (1982)). Methods for Improving the Recovery of Troponin I in Membranes

The adsorption of troponin I and T to surfaces and to various proteins is known to occur and this phenomenon can lower the measured troponin concentration. In particular, when immunoassays are performed in devices or instruments which have a large surface area, for example, when membranes are incorporated into the assay process, the surface area which is exposed to the sample can lower the recovery of troponin. Membranes made up of nylon or compositions of glass fibers having sizes of between 2 mm×2 mm×1 mm and 40 mm×40 mm×5 mm can influence the recovery of troponin I and T when coupled with the assay process.

The troponin I and T molecules have a high degree of basic amino acids. At physiological pH, the basic amino acids are largely positively charged and these charged groups contribute to the adsorptive behavior of the proteins.

In a preferred embodiment of this invention, various components are added to membranes to improve the recovery of troponin I and T in the immunoassay process. Specifically, peptides or proteins which are also strongly basic are added to membranes or surfaces of devices involved in the assay process prior to addition of the sample or reaction mixture. Preferred compounds for this use include peptides, proteins and polymers with pI values greater than about 8. Included in this group are lysozyme, cytochromes, protamine, polylysine, polyvinyl amine, melittin and mastoparan. Concentrations of blocking reagents which are added to surfaces or membranes range from about 0.01 mg/ml to 100 mg/ml and typically about 0.1 mg/ml to 10 mg/ml.

EXPERIMENTAL SECTION

EXAMPLE 1

Preparation of Reagents for Troponin ELISA

Immunoassays

Preparation of Anti-Troponin Antibody Alkaline Phosphatase Conjugates

Alkaline phosphatase (Calzyme, San Luis Obispo, Calif.) in 50 mM potassium phosphate, 10 mM potassium borate, 150 mM sodium chloride, pH 7.0, at 10 mg/ml was derivatized with SMCC (succinimidyl 4-[N-maleimidomethyl] cyclohexane-1-carboxylate, Pierce Chemical Co., Rockford, Ill.) to achieve a 15/1 molar ratio of SMCC/enzyme. The derivatization was performed at room temperature for 90 min and subsequently chromatographed on a GH-25 column (Amicon Corp., Beverly, Mass.) equilibrated in 50 mM potassium phosphate, 10 mM potassium borate, 150 mM sodium chloride, pH 7.0.

The anti-troponin antibodies in 50 mM potassium phosphate, 10 mM potassium borate, 150 mM sodium chloride, pH 7.0, at 10 mg/ml were derivatized with SPDP (N-succinimidyl-3-[2-pyridylthio]propionate, Pierce Chemical Co.) to achieve a 10/1 molar ratio of SPDP/antibody. The antibody was diluted to 2 mg/ml and Dithiothreitol and taurine were added to the solution at final concentrations of 1 mM and 20 mM, respectively, and the solution was subsequently incubated at room temperature for 30 min. The antibody-SPDP was chromatographed on a

GH-25 column (Amicon Corp.) equilibrated in 50 mM potassium phosphate, 10 mM potassium borate, 150 mM sodium chloride, 0.1 mM ethylenediamine tetraacetic acid, pH 7.0.

The SMCC-alkaline phosphatase (in 50 mM potassium phosphate, 10 mM potassium borate, 150 mM sodium chloride, 5 mM magnesium chloride, pH 7.0) and the thiol-antibody (in 50 mM potassium phosphate, 10 mM potassium borate, 150 mM sodium chloride, 0.1 mM ethylenediamine tetraacetic acid pH 7.0), both diluted to 1 mg/ml, were rapidly added to each other with mixing in equimolar amounts. The solution was incubated at room temperature for 3 hours, after which N-ethyl maleimide was added to a final concentration of 2 mM.

Preparation of Biotinylated Troponin Antibodies

Biotin-XX, succinimidyl ester (6-((6(biotinoyl)amino)hexanoyl)amino)hexanoic acid, succinimidyl ester, Molecular Probes, Eugene, Oreg.) at 40 mM in dimethylformamide was added slowly with mixing to an antibody solution at 2 mg/ml in 50 mM potassium borate, 150 mM sodium chloride, pH 8.2, (BBS) to achieve a final molar ratio of 20/1 biotin-XX/antibody. The solution was incubated at room temperature for 2 h, after which the solution was dialyzed at 40° C. for at least 12 h.

Preparation of Avidin-HS Magnetic Latex

Add, while vortexing, 1 ml Estapor Paramagnetic latex particles (0.94 μ , Bangs Laboratories, Carmel, Ind., at 10% solids, washed 4 times with deionized water) in water to 9 ml of 0.55 mg/ml avidin-HS (Scripps Laboratories, San Diego, Calif.) in 50 mM Tris hydrochloride, 150 mM sodium chloride, pH 7.5. Incubate the latex solution at 45° C. for 2 h. Wash the latex 3 times, each with 10 ml BBS, and resuspended in 10 ml BBS.

EXAMPLE 2

Immunoassay of Human cardiac Troponin I and Troponin T

Two immunoassay protocols are described. They were used to detect Troponin I and Troponin T, either present in human serum and plasma or in solutions containing purified proteins.

Protocol A

The sample containing troponin I or troponin T was diluted to 1-10 ng/ml troponin I or troponin T in an assay buffer (hereafter called assay buffer) containing 10 mM 3-(N-morpholino)propane sulfonic acid, 650 mM sodium chloride, 1 mM magnesium chloride, 0.1 mM zinc chloride, 1 mg/ml polyvinyl alcohol (10,000 m.w.), 10 mg/ml bovine serum albumin, 1 mg/ml sodium azide, pH 7.0. To 25 μ l of diluted sample in a microtiter plate well was added 50 μ l of assay buffer containing 2.5 μ g/ml anti-troponin I or troponin T antibody conjugates (Example 1) and 2.5 μ g/ml biotinylated anti-troponin I or troponin T polyclonal antibody (Example 1) to form a reaction mixture. After a 30 minute incubation of the reaction mixture at room temperature, 25 μ l of avidin-HS coated magnetic latex (Example 1; 0.5% latex in assay buffer) was added to the microtiter plate well, followed by a 5 minute incubation at room temperature. The magnetic latex was pelleted using a microtiter plate magnet (Perceptive Diagnostics, Cambridge, Mass.) and washed twice in BBS-Tween (20 mM borate, 150 mM sodium chloride, 0.1 mg/ml sodium azide, 0.02% Polyoxyethylene-20-Sorbitan Monolaurate (Tween-20), pH 8.2) and once in TES (40 mM Tris, 150 mM sodium chloride, pH 7.5) The pellet was resuspended in ELISA amplification reagents (Gibco BRL, Gaithersburg, Md.) according to the manufacturer's instructions. After the amplification was complete,

the magnetic latex was pelleted and 80 μ l of the colored supernatant was transferred to a fresh microtiter plate. The absorbance at 490 nm (OD_{490}) was measured using a microtiter plate reader (Molecular Devices, Palo Alto, Calif.).

Protocol B

The sample containing troponin I or troponin T was diluted into assay buffer as described in protocol A. To 80 μ l of diluted sample in a microtiter plate well was added 40 μ l of assay buffer also containing 30 μ g/ml anti-troponin I or troponin T monoclonal antibody and 7.5 μ g/ml biotinylated anti-troponin I or troponin T polyclonal antibody (Example 1) that was complimentary to the monoclonal antibody to form the reaction mixture. Aliquots (25 μ l) were removed at various times (2 minutes to 24 hours) and were added to microtiter plate wells containing 25 μ l of avidin-HS coated magnetic latex (0.5% latex solids in assay buffer), followed by a 5 minute incubation. The magnetic latex was pelleted and washed once in BBS-Tween and once in assay buffer. The pellet was resuspended in 25 μ l of assay buffer also containing 5 μ g/ml of goat anti mouse kappa antibody conjugated to alkaline phosphatase (Southern BioTechnology Associates, Inc., Birmingham, Ala.) followed by a 15 minute incubation. The magnetic latex was pelleted and the remainder of the assay was performed as indicated in Protocol A.

EXAMPLE 3

Selection of Anti-Troponin I Antibodies that Bind the Oxidized, the Reduced or Both the Oxidized and Reduced Forms of Human Cardiac Troponin I

Anti-troponin I antibody conjugates (Example 1) and complimentary biotinylated troponin I polyclonal antibodies (Example 1) were tested for recognition of intramolecularly oxidized or reduced troponin I. The anti troponin I monoclonal antibodies tested were: clone 2D5 and clone 1A12 (BiosPacific, Emeryville, Calif.), clone 110 and 111 (Research Diagnostics, Inc., Flanders, N.J.) and clone TRI-7 F81 (DAKO Corporation, Carpinteria, Calif.). The biotinylated anti-troponin I antibodies tested were affinity-purified goat polyclonals, specified as peptide 1, peptide 2, peptide 3 or peptide 4 specific (BiosPacific, Emeryville, Calif.). Human cardiac troponin I (P. Cummins, University of Birmingham, Birmingham, UK) was air oxidized at 1.0 μ g/ml as described in example 4 to form the intramolecular disulfide. The oxidized troponin I was diluted to 1–10 ng/ml in assay buffer either without (oxidized sample) or with (reduced sample) dithiothreitol (DTT) at a final concentration of 3mM, followed by a 3 hour incubation at room temperature to allow reduction of the disulfide by DTT. The oxidized and reduced samples were assayed without further dilution using Protocol A of Example 2. The results are shown in Table 1 and are expressed in terms of a ratio of the slope for oxidized troponin I (TnI) divided by the slope for reduced troponin I. The assay slope increases with increasing recognition of troponin I by the antibody pair.

The data show that antibodies can be selected that either preferentially bind oxidized or reduced troponin I or bind oxidized and reduced troponin I approximately the same. Selection of the antibodies without regard to the oxidation-reduction state of troponin I can lead to a substantial error in the quantification of the troponin I concentration.

TABLE 1

	Anti troponin I monoclonal antibody	Anti troponin I polyclonal antibody	Ratio of assay slopes (oxidized TnI/reduced TnI)
5	Clone 2D5	peptide 1 specific	8.3
10	Clone 2D5	peptide 3 specific	10
	Clone 1A12	peptide 1 specific	0.6
	Clone 1A12	peptide 3 specific	1.3
15	Clone 1A12	peptide 4 specific	1.2
	Clone TRI-7 F81	peptide 1 specific	0.5
	Clone TRI-7 F81	peptide 2 specific	0.5
20	Clone TRI-7 F81	peptide 3 specific	0.5
	Clone TRI-7 F81	peptide 4 specific	0.5
	Clone 110	peptide 4 specific	0.8
25	Clone 111	peptide 3 specific	1.0
	Clone 111	peptide 4 specific	0.4

EXAMPLE 4

Oxidation-Reduction of Purified Human Cardiac Troponin I
The kinetics of intramolecular oxidation and reduction of purified troponin I (P. Cummins, University of Birmingham, UK) was measured with an immunoassay (Protocol A, Example 2) using a clone 2D5 antibody conjugate (Example 1) and biotinylated goat anti troponin I peptide 1 polyclonal antibody (Example 1). This antibody pair binds strongly to oxidized troponin I and weakly to reduced troponin I as described in Example 3. The results of the assay are expressed in terms of an assay slope [OD_{490} per ng/ml total (oxidized+reduced) troponin I] in the linear range of the assay. The assay slope increases with the fraction of oxidized troponin I.

Air Oxidation of Reduced Troponin I

The rate of air oxidation of troponin I at two troponin I concentrations was measured. Reduced troponin I at 0.27 mg/ml in a buffer containing 20 mM Tris-HCl, 0.5M sodium chloride, 60 mM 2-mercaptoethanol, pH 7.5, was diluted to either 1300 ng/ml or 10 ng/ml in assay buffer containing either no or 25 mM 2-mercaptoethanol. The solutions were incubated at room temperature. Aliquots were taken after various incubation times, as indicated in FIG. 1a, diluted to 4 and 8 ng/ml troponin I in assay buffer, and assayed immediately. The results are shown in FIG. 1a, wherein the error bars represent 1 standard deviation (SD).

Peroxide Oxidation of Reduced Troponin I

Peroxide (Fisher, unstabilized) was added to final concentration of 2 mM or 20 mM to an aliquot of the 1300 ng/ml solution of reduced troponin I (see this example, air oxidation) immediately after the troponin I was diluted from the 0.27 mg/ml stock solution. Aliquots were taken after various incubation times, as indicated in FIG. 1b, treated with catalase (Calbiochem, La Jolla, Calif.; 0.01 mg/ml final concentration for 5 minutes) to remove the peroxide, diluted

to 4 and 8 ng/ml troponin I, and assayed immediately. The results are shown in FIG. 1b, wherein the error bars represent 1 SD.

DTT Reduction of Oxidized Troponin I

Troponin I that was incubated (air oxidized) at 1000 ng/ml in assay buffer for 15 hours at room temperature was diluted to 4 and 8 ng/ml in assay buffer. DTT was added to a final concentration of 0, 1.5 and 3.0 mM followed by incubation at room temperature for the times indicated in FIG. 2. The aliquots were then assayed for troponin I. After steady state was reached (approximately 6 hours), aliquots (100 μ l) from the three DTT concentration samples were reoxidized with 20 mM peroxide for 15 minutes, treated with catalase for 5 minutes and assayed. The results are shown in FIG. 2, wherein the error bars represent 1 SD.

The data show that the results of an immunoassay can vary over time if the oxidation-reduction state of the troponin I is allowed to change. The oxidation-reduction state of troponin I, and thus the immunoassay results, can be reversibly changed and greatly stabilized over time by the use of oxidants and reductants.

EXAMPLE 5

Alkylation of Reduced Troponin I

Troponin I was rapidly alkylated using various alkylating reagents. The stock reduced troponin I (University of Birmingham) was at 0.27 mg/ml in 20 mM Tris hydrochloride, 0.5M sodium chloride, 50 mM 2-mercaptoethanol. Three alkylation reactions (#1-3) were performed and a control was prepared (#4):

1. 20 μ l of stock troponin I was added to 20 μ l 0.5M potassium borate, 0.2 mM ethylenediamine tetraacetic acid, pH 8.0 and subsequently, 10 μ l 398 mM iodoacetamide was added.
2. 20 μ l of stock troponin I was added to 20 μ l 0.5M potassium borate, 0.2 mM ethylenediamine tetraacetic acid, pH 8.0 and subsequently, 10 μ l 398 mM iodoacetic acid was added.
3. 20 μ l of stock troponin I was added to 20 μ l 0.5M potassium borate, 0.2 mM ethylenediamine tetraacetic acid, pH 8.0 and subsequently, 12.5 μ l 319 mM N-ethylmaleimide was added.
4. 20 μ l of stock troponin I was added to 20 μ l 0.5M potassium borate, 0.2 mM ethylenediamine tetraacetic acid, pH 8.0 and subsequently, 10 μ l 0.5M potassium borate, 0.2 mM ethylenediamine tetraacetic acid, pH 8.0 was added.

The reactions were incubated at room temperature for 1 h 25 min. During this incubation, the stock troponin I was kept on ice. Aliquots (24 μ l) of each solution (1-4) were added to 0.9 μ l of 2.9M mercaptoethanol and were incubated at room temperature for 15 min, after which the samples were frozen in liquid nitrogen. The remaining aliquots of each solution (1-4) were also frozen in liquid nitrogen with no further treatment.

EXAMPLE 6

Immunoassay of Alkylated Troponin I

Freshly thawed Troponin I alkylated (Example 5) with N-ethyl maleimide (NEM), iodoacetic acid (IHAC), iodoacetamide (IAM), or no alkylation reagent (control sample, Example 5) was diluted to 1-10 ng/ml in assay buffer. A freshly thawed aliquot of the reduced stock troponin I (Example 5) was diluted (standard sample) into assay buffer containing either 0 or 3 mM DTT. Aliquots (25 μ l) of all dilutions were taken after either a 0.5 hour or 5.5 hour incubation at room temperature and assayed (Protocol A,

Example 2) using a clone 2D5 anti troponin I monoclonal antibody conjugate and biotinylated goat anti troponin I peptide 1 specific polyclonal antibody (Example 1). This antibody pair binds strongly to oxidized troponin I and weakly to reduced troponin I (Examples 3 and 4). Troponin I will remain substantially reduced during a 0.5 hour incubation but will be almost completely oxidized (by air) after a 5.5 hour incubation unless DTT is present to stabilize the reduced form (see Example 4). The results are shown in Table 2 and are given in terms of assay slope (OD₄₉₀ per ng/ml total troponin I) in the linear range. A larger assay slope indicates a stronger binding interaction between the antibodies and the troponin I.

The data show that the antibody pair binds alkylated troponin I similarly to reduced troponin I, that is, weakly in comparison with oxidized troponin I. Furthermore, alkylation stabilizes the immunoassay result with respect to time, similarly to the effect observed by the use of oxidants or reductants to stabilize the oxidation-reduction state of troponin I (Example 4). The lower and more stable assay slope of the control sample as compared with the standard sample is explained by the presence of mixed disulfides formed between the two cysteine residues of the control sample troponin I and 2-mercaptoethanol during the room temperature incubation of the control sample at pH 8 (see Example 5).

TABLE 2

Sample	Assay slope (0.5 hour incubation)	Assay slope (5.5 hour incubation)	Ratio of assay slopes (5.5 hour/0.5 hour)
reduced TnI standard (+DTT)	0.030	0.030	1.0
reduced TnI standard (-DTT)	0.058	0.21	3.6
TnI Control	0.037	0.078	2.1
TnI alkylated with NEM	0.023	0.026	1.1
TnI alkylated with IAM	0.013	0.013	1.0
TnI alkylated with IHAC	≤ 0.01	≤ 0.01	

EXAMPLE 7

Effect of Peroxide on Immunoassay of Cardiac Troponin I from Patients with Confirmed Myocardial Infarction

Frozen Human serum or plasma, drawn in heparin tubes from patients with confirmed myocardial infarction, was obtained from local hospitals. The serum or plasma was thawed at room temperature and immediately split into two aliquots. One aliquot was oxidized at room temperature by the addition of peroxide at a final concentration of 20 mM. The second aliquot was untreated. The oxidation reaction was stopped after 20 minutes by the addition of catalase at a final concentration of 0.01 mg/ml. Ten minutes after the catalase was added both the oxidized and the untreated aliquots were diluted serially by factors of four in assay buffer and assayed immediately for cardiac troponin I using the 2D5 anti troponin I conjugate and biotinylated anti troponin I peptide 3 specific antibodies (Example 2, Protocol A). This complimentary antibody pair binds oxidized troponin I strongly and reduced troponin I weakly (Example 3). Air oxidized (example 4) purified troponin I (P. Cummins,

University of Birmingham), diluted to 2, 4, and 8 ng/ml in assay buffer, was assayed with the same antibody reagents to construct a standard curve. The concentration of troponin I in the neat oxidized or untreated serum or plasma sample (Table 3) was calculated from this standard curve using the OD₄₉₀ measurements that fell within the linear range of the assay.

The data show that oxidation of serum or plasma samples from patients with confirmed myocardial infarction can have a substantial effect on the concentration of cardiac troponin I determined by immunoassay. Immunoassay of troponin I in serum or plasma without regard to the oxidation state of the troponin I could lead to a serious underestimation of the troponin I concentration and result in the non-diagnosis of a myocardial infarction.

TABLE 3

Sample	Time between sample collection and freezing (hours)	Troponin I concentration by assay of untreated sample (ng/ml)	Troponin I concentration by assay of Peroxide oxidized sample (ng/ml)
1 Plasma	2	6.3	9.4
2 Plasma	6.5	0.8	1.0
3 Serum	9.3	6.6	8.3
4 Serum	6.5	31.9	46.5
5 Plasma	6.5	31.7	49.7
6 Plasma	9.5	0.6	1.0
7 Serum	11.5	0.4	0.4
8 Plasma	5.0	4.5	5.4
9 Serum	10.5	1.6	2.3
10 Plasma or Serum	unknown	13.6	13.2

EXAMPLE 8

Effect of Peroxide on Immunoassay of Cardiac Troponin I in Human Plasma after two Freeze/Thaw Cycles

Plasma sample number 5 (Table 3, Example 7) was stored untreated on ice for three hours after it was initially thawed and then refrozen and stored at -70° C. for several days. The plasma was thawed at room temperature and split into two aliquots; one was oxidized with peroxide and the other was left untreated as described in Example 7. The concentration of troponin I in the oxidized and untreated aliquots was determined immediately by the immunoassay described in example 7 and was found to be 53.9 ng/ml in the untreated aliquot and 56.4 ng/ml in the oxidized aliquot.

The data show that oxidation of the plasma after the second thaw did not have a substantial effect on the concentration of cardiac troponin I determined by immunoassay.

EXAMPLE 9

Immunoassay of Cardiac Troponin I in Oxidized and Reduced Plasma From a Patient With Myocardial Infarction

Frozen Human plasma drawn in heparin tubes from a patient with a confirmed myocardial infarction was obtained from a local hospital. The plasma was thawed at room temperature and immediately split into two aliquots. One aliquot was oxidized with peroxide as described in Example 7. The other aliquot was reduced by addition of DTT to a final concentration of 10 mM, followed by a 3 hour incu-

bation at room temperature. The oxidized aliquot was then diluted serially by factors of 2 into assay buffer and the reduced aliquot was diluted serially by factors of 2 into assay buffer containing 3 mM DTT. The diluted aliquots were assayed for troponin I immediately (Protocol A, Example 2) either with the complementary antibody pair clone 2D5 anti troponin I conjugate and biotinylated anti troponin I peptide 3 polyclonal antibody or with the complementary antibody pair clone TRI-7 F81 anti troponin I conjugate and biotinylated anti troponin I peptide 3 polyclonal antibody. Purified troponin I (University of Birmingham) which was air oxidized (example 4) was diluted to 2, 4, and 8 ng/ml in assay buffer and was used to construct the standard curve from which the concentration of troponin I in the neat oxidized or reduced plasma sample was determined. The results are shown in Table 4.

The data show that chemical oxidation and reduction of cardiac troponin I in the plasma sample affects the recognition of the tested antibody pairs for the troponin I in a manner similar to that observed for purified troponin I (Example 3).

TABLE 4

Monoclonal antibody conjugate	Assayed troponin I concentration (ng/ml) in oxidized plasma	Assayed troponin I concentration (ng/ml) in reduced plasma	Ratio of troponin I concentrations (oxidized plasma/reduced plasma)
Clone 2D5	82	<1	>82
Clone TRI-7 F81	52.8	77.5	0.68

EXAMPLE 10

Selection of Anti Troponin I Antibodies That Are Either Sensitive or Insensitive To The Binding of Troponin C to Troponin I

Monoclonal anti troponin I conjugates and complimentary biotinylated anti troponin I polyclonal antibodies (Example 1) were tested for their recognition of free troponin I and troponin II bound to troponin C in a binary complex. Four types of troponin I samples were prepared at room temperature and assayed for troponin I; they are: oxidized troponin I with and without added troponin C and reduced troponin I with and without added troponin C. Oxidized (by air, see Example 4) Human cardiac troponin I (P. Cummins, University of Birmingham) was diluted to 2, 4, and 8 ng/ml in assay buffer containing 2 mM calcium chloride. One aliquot of each concentration of troponin I was either untreated or reduced by the addition of DTT to a final concentration of 3 mM from a 30 mM DTT stock solution in assay buffer to form a reduction reaction. Three hours after the reduction reaction was started, each oxidized and reduced troponin I aliquot was split into two aliquots; to one aliquot was added human cardiac troponin C (Bio-Tech International Inc., Seattle, Wash.) to a final concentration of 0.1 mg/ml from a 1 mg/ml stock solution in 20 mM potassium phosphate, 4 mM potassium borate, 150 mM sodium chloride, pH 7.0 to form a binding reaction mixture, and to the other aliquot was added the same volume of the above buffer without troponin C. One hour after the troponin C was added, all the aliquots were assayed for troponin I (Protocol A, Example 2) using the antibody pairs listed in Table 5. The results in Table 5 are expressed as a fractional assay response which was determined by dividing the assay

slope in the presence of troponin C by the assay slope in the absence of troponin C.

The results in Table 5 show that some antibody pairs recognize free troponin I and troponin I bound to troponin C equally well, while other antibody pairs recognize only free troponin I. An immunoassay with antibodies that do not recognize troponin I bound to troponin C will underestimate the total troponin I concentration when some of the troponin I is present as the troponin I/C binary complex.

TABLE 5

Anti troponin I antibody conjugate	Biotinylated anti troponin I polyclonal antibody	Fractional assay response	
		Oxidized troponin I	Reduced troponin I
Clone 2D5	Peptide 3 specific	0.81	0.60
Clone 111	Peptide 1 specific	0.83	Not determined
Clone 111	Peptide 3 specific	0.47	0.52
Clone 111	Peptide 4 specific	0.59	0.19
Clone 110	Peptide 4 specific	0.96	0.48
Clone 1A12	Peptide 1 specific	<0.05	<0.05
Clone 1A12	Peptide 3 specific	<0.05	<0.05
Clone 1A12	Peptide 4 specific	<0.05	<0.05
Clone TR7 F81	Peptide 1 specific	0.74	0.79
Clone TR7 F81	Peptide 2 specific	0.92	1.04
Clone TR7 F81	Peptide 3 specific	0.94	0.97
Clone TR7 F81	Peptide 4 specific	0.70	0.79

EXAMPLE 11

Effect of Troponin T, EDTA and Melittin on a Troponin I Immunoassay With Troponin C Present in large Excess Over Troponin I

Ethylenediamine tetraacetic acid (EDTA) lowers the binding affinity of troponin I for troponin T and troponin C by chelating calcium and magnesium ions. Melittin lowers the affinity of troponin I for troponin C by binding to troponin C. The effectiveness of EDTA and Melittin (hereafter referred to as binding inhibitors) in breaking up the binary complex of troponin I and troponin C in the presence and absence of troponin T was assessed. Oxidized Human cardiac troponin I (P. Cummins, University of Birmingham) at 1.0 $\mu\text{g/ml}$ in assay buffer containing 2 mM calcium chloride was reduced with dithiothreitol at a final concentration of 3 mM for three hours at room temperature. The reduced troponin I was diluted to 2 and 4 ng/ml in assay buffer containing 2 mM calcium chloride and 3 mM dithiothreitol. Each concentration was split into four aliquots to which were added human cardiac troponin C (Bio-tech International, Inc.) to final concentrations of 0, 0.1, 1.0, and 10.0 $\mu\text{g/ml}$ and from 100-fold excess stock solutions in 20 mM potassium phosphate, 4 mM potassium borate, 150 mM sodium chloride, pH 7.0. Each of the resulting aliquots were further split into two aliquots to which was added human troponin T (Scripps Labs) to a final concentration of either

0.0 or 0.1 $\mu\text{g/ml}$ from a 100-fold excess stock solution in deionized water. The aliquots were incubated at room temperature for one hour after the addition of troponin T, then assayed for troponin I (Protocol B, Example 2). The antibody solution added to the microtiter plate wells contained 30 $\mu\text{g/ml}$ clone 1A12 anti troponin I and 7.5 $\mu\text{g/ml}$ biotinylated anti troponin I peptide 4 specific antibodies (example 1) either without or with binding inhibitors (30 mM EDTA and 0.15 mM Melittin (Sigma Chemical, Co., St. Louis, Mo.)). Aliquots of the reaction mixtures formed by the addition of antibodies to the troponin I samples were removed after 0.5 h and were further treated as described in Protocol B, Example 2. The assay results for the samples containing no troponin C or T and no binding inhibitors were used to construct a standard dose-response curve. The fractional assay response (shown in FIG. 3) for samples containing inhibitors and troponin components was determined by dividing the assay slope for each sample by the slope of the standard curve.

The data show that in the presence of troponin C, the troponin I concentration is largely underestimated. The binding inhibitors almost completely reverse the effect of troponin C. The presence of troponin T in the absence of troponin C has no effect on the troponin I immunoassay. In the presence of troponin C and T, the concentration of troponin I is dramatically reduced. The binding inhibitors appear to be less effective at opening up or partially unraveling the troponin complex when the complex is ternary than when the complex is binary.

EXAMPLE 12

Effect of Binding Inhibitors on an Immunoassay of Troponin I in the Presence of Troponin C or Troponin C and T

Solutions containing 1.0 $\mu\text{g/ml}$ purified human cardiac troponin I (reduced by DTT, Example 4) and either 1.2 $\mu\text{g/ml}$ human cardiac troponin C (Bio-tech International, Inc.) or 1.2 $\mu\text{g/ml}$ troponin C and 3.1 $\mu\text{g/ml}$ human cardiac troponin T (Scripps Labs) were incubated for 2 hours at room temperature in assay buffer containing 3 mM DTT. Troponin C was added to troponin I prior to addition of troponin T. The troponin solutions were diluted to 2, 4 and 8 ng/ml in terms of troponin I concentration in assay buffer containing 2 mM calcium chloride and 0.5 mM DTT and assayed immediately with and without binding inhibitors as described in example 11. Aliquots of the reaction mixtures of antibodies and troponin components were removed 0.5 h and 2.2 h after the antibodies were added. These aliquots were further treated as described in Protocol B, Example 2. Reduced troponin I without added troponin C and T and without binding inhibitors was assayed to produce a standard curve. The results are expressed in Table 6 as a fractional assay response which was determined by dividing the assay slope for each sample by the slope of the standard curve.

The data show that troponin T and C present in approximately a two fold molar excess above the troponin I concentration substantially lowers the amount of troponin I measured in the immunoassay. Troponin C alone has a smaller effect on the measured troponin I concentration at the antibody concentrations used in this assay. The binding inhibitors partially reverse the effect of troponin C and T at 0.5 h incubation and completely at 2.2 h.

TABLE 6

Time after antibodies added (hours)	Fractional assay response			
	With troponin C		With troponin C and T	
0.5	0.88	0.94	0.55	0.75
2.2	1.15	1.02	0.44	1.02

EXAMPLE 13

Assay of Purified Human Cardiac Ternary Troponin Complex for Troponin I

Purified human cardiac ternary troponin complex (Bio-tech International, Inc., 3 mg/ml stock solution in a buffer containing 20 mM sodium citrate, pH 6) was diluted to concentrations ranging from 1 to 15 ng/ml total troponin I in assay buffer containing 2 mM calcium chloride but without the 0.1 mM zinc chloride or the 1 mM magnesium chloride (hereafter called metal-free assay buffer). The diluted solutions were assayed for troponin I with and without inhibitors present as described in Example 11. Aliquots of the reaction mixtures of the antibodies with the troponin complex were taken at the times indicated in FIG. 4 and were further treated as described in Protocol B, Example 2. Purified troponin I (Bio-tech International, Inc.) was assayed and used to construct a standard curve. The fractional assay response shown in FIG. 4 was determined by dividing the assay slope for the complex (OD_{490} as a function of total troponin I concentration) by the slope of the standard curve.

The results show that troponin I that is bound in the ternary complex with troponin C and T is not recognized very well by the antibodies at the antibody and troponin I concentrations used in this assay, even after a very long incubation time. In particular, the 30 minute time point had no detectable troponin I, with or without inhibitors. The binding inhibitors slowly open up or partially unravel the complex and thus slowly increase the fraction of troponin I recognized by the antibodies.

EXAMPLE 14

Effect of Binding Inhibitors on a Troponin I Immunoassay of Plasma and Serum from Patients with Confirmed Myocardial Infarction

Frozen serum and plasma drawn in EDTA and Heparin tubes from patients with a confirmed myocardial infarction were obtained from local hospitals and thawed at room temperature. Calcium chloride was added to a final concentration of 6 mM (to bind all the EDTA) and the resultant solution was incubated for two hours at room temperature and then overnight at 4° C. The incubated samples were diluted by factors of two to a maximum dilution factor of 256 in metal-free assay buffer. The diluted samples were immediately assayed for troponin I with and without binding inhibitors as described in Example 11, with the exception that the polyclonal antibody was goat anti-troponin I peptide 3 specific. Aliquots of the reaction mixture formed by the addition of antibodies to the diluted samples were taken at the times indicated in FIGS. 5a-f and further treated as described in Example 2, Protocol B. Oxidized troponin I (University of Birmingham) at 2.4 and 8 ng/ml in metal-free assay buffer was assayed to produce a standard curve. The OD_{490} values measured for the diluted serum or plasma

samples were plotted as a function of the inverse of the dilution factor. The slope in the linear region of the resultant curve (typically at $OD_{490} < 2$, which corresponds to a troponin I concentration of less than 8 ng/ml) was divided by the slope of the standard curve to obtain the concentration of troponin I in the neat serum or plasma sample shown in FIGS. 5a-f. Each FIG. 5a through 5f reflects immunoassays on serum or plasma from different patients.

The data show that the measured concentration of troponin I in all of the serum and plasma samples tested was increased by the addition of binding inhibitors. Importantly, the time profile of the measured concentration of troponin I was in some cases biphasic (FIGS. 5a-c and 5e). The fast phase was complete within the first assay time of 0.5 h. The slow phase continued to rise for 6-24 hours depending on the sample. A slow phase was also observed for the purified ternary troponin complex (Example 13). The slow phase observed for the diluted serum and plasma samples, may, therefore, be associated with the opening up or partial unraveling of a ternary complex by the inhibitors and antibodies. The fast phase indicates a bound complex of troponin I that is more easily opened up or partially unraveled than the ternary complex, since the fast phase is absent for the purified ternary complex (Example 13). Thus, the fast phase could be associated with the opening up or partial unraveling of binary complexes of troponin I.

EXAMPLE 15

Immunoassays That are Sensitive To Free Troponin I, Troponin I Bound in a Ternary Complex, and Both Free and Bound Troponin I

Three sets of antibodies were evaluated for their ability to recognize free troponin I and troponin I bound in the ternary complex. Three antibody stock solutions described below as #1-3 were prepared in metal-free assay buffer either with or without binding inhibitors (30 mM EDTA and 0.15 mM Melittin) and the following antibodies:

1. 30 µg/ml 1A12 anti troponin I and 7.5 µg/ml biotinylated anti troponin I peptide 4 specific;
2. 30 µg/ml 1A12 anti troponin I, 30 µg/ml 9B1 anti troponin T monoclonal (Biospecific), 5 µg/ml each of biotinylated anti troponin I peptide 1, 2, 3 and 4 specific;
3. 30 µg/ml 9B1 anti troponin T and 5 µg/ml each of biotinylated anti troponin I peptide 1, 2, 3 and 4 specific.

Human cardiac troponin ternary complex (Bio-tech International, Inc.) was diluted to 1-15 ng/ml troponin I equivalents in metal free assay buffer and purified troponin I (Bio-tech International, Inc.) was diluted to 2.4 and 8 ng/ml in the same buffer. The dilutions of the troponin complex and troponin I were assayed immediately using Protocol B, Example 2 with antibody solutions #1-3. Aliquots of the reaction mixtures formed by the addition of the antibodies to the troponin complex and troponin I samples were taken at 0.5 h and 2.5 h after the antibodies were added and were further treated as described in Example 2, Protocol B. The results are shown in Table 7 and are expressed in terms of an assay slope with units of OD_{490} per ng/ml total troponin I in the linear range of the assay. A higher slope indicates better binding of the antibodies to the troponin components. Antibody solution #2 was tested and found to be negative for cross reactivity with purified human cardiac troponin T (Scripps Labs) at 1-6 ng/ml using the assay protocol described herein (data not shown).

The data show that the immunoassay using the antibodies in solution #1 recognizes free troponin I but not troponin I in the ternary complex. The immunoassays using the antibodies in solution #2 recognizes free troponin I and troponin I in the ternary complex almost equally well. Thus, antibody solution #2 is superior to solution #1 for the assay of total troponin I when a fraction of the troponin I is present as the ternary complex. The immunoassay using the antibodies in solution #3 recognizes the ternary complex well but recognizes free troponin I poorly. This poor recognition of free troponin I causes the assay slope for antibody solution #3 to decrease over time in the presence of binding inhibitors. By using all three antibody solutions in immunoassays, one can estimate independently the concentrations of free troponin I (solution #1), total troponin I (solution #2) and bound troponin I (solution #3).

TABLE 7

Antibody solution #	Time after antibodies added	Assay Slope (OD ₄₉₀ per ng/ml total troponin I)		
		Free Troponin I	Troponin ternary complex (without binding inhibitors)	Troponin ternary complex (with binding inhibitors)
#1	0.5	0.14	<0.006	<0.006
#1	2.5	0.14	0.003	0.021
#2	0.5	0.12	0.18	0.14
#2	2.5	0.22	0.17	0.12
#3	0.5	0.008	0.18	0.17
#3	2.5	0.008	0.17	0.08

EXAMPLE 16

Estimation of Free Troponin I, Bound Troponin I and Total Troponin I in Plasma From a Patient with a Myocardial Infarction

The three antibody solutions described in Example 15 were used in immunoassays to measure the troponin I concentration in plasma from a patient with a confirmed myocardial infarction. The frozen plasma was treated and diluted into metal-free assay buffer as described in Example 14. The diluted plasma was assayed for troponin I immediately using Protocol B, Example 2 with antibody solutions #1-3 (Example 15). Aliquots were taken 0.5 h and 2.5 h after the antibodies were added to the samples to form the reaction mixtures and were further treated as described in Example 2, Protocol B. Either free troponin I at 1-4 ng/ml or the troponin complex (Bio-tech International, Inc.) at 1-8 ng/ml troponin I in metal-free assay buffer was assayed and used to construct a standard curve of OD₄₉₀ as a function of total troponin I concentration for each antibody solution. The concentration of troponin I in the neat plasma sample, as shown in Table 8, measured by each antibody solution was determined using either the standard curve for free troponin I or for troponin I in the ternary complex in the absence of binding inhibitors as indicated in Table 8 and as described in Example 14. The ratio determined by dividing the assayed troponin I concentration with binding inhibitors by the concentration without inhibitors is also shown in Table 8.

The data show that the concentration of troponin I determined by immunoassay using antibody solution #1 is much more sensitive to the presence of binding inhibitors and thereby to the opening up or partial unraveling of the

troponin complex than that determined using solution #2 or #3. The data of Example 15 suggest that antibody solution #1 measures mainly free troponin I, solution #2 measures both free troponin I and troponin I bound in the ternary complex and solution #3 measures mainly troponin I bound in the ternary complex. Thus, the conclusions from Example 15 taken together with the data in Table 8, indicate that substantial amounts of both free and bound troponin I are present in the diluted plasma sample. Among the three antibody solutions, used in the immunoassays, solution #2 gives the largest assayed troponin I concentration, as expected, because the immunoassay using antibody solution #2 measures both free and bound troponin I. Thus, antibody solution #2 appears to provide the most sensitive measure of troponin I in the plasma sample. Antibody solution #2 standardized with the purified troponin complex gave the most stable assay with respect to inhibitor addition and assay incubation time. The decrease of troponin I concentration at 2.5 h measured with antibody solution #2 standardized with free Troponin I is due to an increase at 2.5 h of the slope of the standard curve.

TABLE 8

Antibody solution	Time after antibodies added, (hours)	Standard used to determine troponin I concentration	Troponin I concentration in plasma, (ng/ml)		Ratio
			Without binding inhibitors	With binding inhibitors	
#1	0.5	Troponin I	91	186	2.0
#1	2.5	Troponin I	83	251	3.0
#3	0.5	Ternary Troponin Complex	134	108	0.8
#3	2.5	Ternary Troponin Complex	87	65	0.74
#2	0.5	Troponin I	280	360	1.3
#2	2.5	Troponin I	172	220	1.3
#2	0.5	Ternary Troponin Complex	229	299	1.3
#2	2.5	Ternary Troponin Complex	223	286	1.3

EXAMPLE 17

Immunoassay of Free Human Cardiac Troponin T and Troponin T in the Human Cardiac Ternary Complex

Two antibody stock solutions (#1 and 2) were prepared as described below in metal-free assay buffer either with or without binding inhibitors (30 mM EDTA and 0.15 mM Melittin) and with the following antibodies:

- 30 µg/ml 1A12 anti troponin I, 30 µg/ml 9B1 anti troponin T and 7.5 µg/ml biotinylated anti troponin T peptide 3 specific (Biospecific).
- 30 µg/ml 9B1 anti troponin T and 7.5 µg/ml biotinylated anti troponin T peptide 3 specific.

Human cardiac troponin ternary complex (Bio-tech International, Inc.) was diluted to 1-20 ng/ml troponin T in metal free assay buffer and purified Human cardiac troponin T (Scripps Labs) was diluted to 1.5, 3.0 and 6.0 ng/ml in the

same buffer. The dilutions of the troponin complex and troponin T were assayed immediately using Protocol B, Example 2 with antibody solutions #1 and #2. Aliquots of the reaction mixtures formed by the addition of the antibodies to the troponin complex and troponin T samples were taken at 0.5 h and 3.0 h after the antibodies were added and were further treated as described in Example 2, Protocol B. The results are shown in Table 9 and are expressed in terms of an assay slope with units of OD₄₉₀ per ng/ml total troponin T in the linear range of the assay. A higher slope indicates better binding of the antibodies to the troponin components.

The data show that the antibodies in solution #1 recognize both free troponin T and troponin T bound in the ternary complex equally well. The antibodies in solution #2 recognize free troponin T well but recognizes troponin T in the ternary complex poorly. Thus, antibody solution #1 is expected to provide the most sensitive and accurate measure of the total concentration of troponin T in Human blood samples in which a substantial amount of the ternary complex is present.

TABLE 9

Antibody solution	Time after antibodies added to troponin sample (hours)	Assay slope (OD ₄₉₀ per ng/ml troponin T)			
		Free troponin T without binding inhibitors	Free troponin T with binding inhibitors	Troponin complex without binding inhibitors	Troponin complex with binding inhibitor
#1	0.5	0.069	0.078	0.061	0.070
#2	0.5	0.078	0.085	0.013	0.012
#2	3.0	>0.08	>0.08	0.018	0.023

EXAMPLE 18

Use of Troponin C to Prevent Non-specific Binding of Troponin I to Filter Membranes

oxidized cardiac troponin I (air oxidized, Example 4, University of Birmingham) at 100 ng/ml final concentration in human serum (Hybritech, Inc., San Diego) was incubated either with or without 100 ug/ml human cardiac troponin C (Bio-tech International, Inc.) at room temperature for 30 minutes. Two filter membranes, a CytoSep filter (Ahlstrom Filtration, Mount Holly Springs, Pa.) and a glass fiber filter (GB-100R, Micro Filtration Systems, Dublin, Calif.) were cut into rectangles measuring 1.5 cm by 3.0 cm and were secured to a transparency film (catalog #pp2500, 3M, Austin, Tex.) by a piece of tape across the filters. The troponin I solutions with and without troponin C (300 µl) were applied slowly to the top of the filters at one end and the solution migrated through the filter to the far end by wicking action. About 15 µl of solution was collected from the far end with a plastic pipet tip. The collected solutions were diluted by factors of 20, 40 and 80 in assay buffer and assayed using Protocol A, Example 2 with TRI-7 F81 anti troponin I conjugate and biotinylated anti troponin I peptide 3 specific antibodies. Aliquots of the troponin I solutions with and without troponin C that had not been passed through the filters were also assayed and used to construct a standard curve from which the concentration of troponin I in the solution that had passed through the membrane was determined. The calculated concentration was divided by 100 ng/ml to obtain the fraction of recovered troponin I shown in Table 10. The experimental errors given in Table 10 represent one standard deviation.

The data show that the presence of troponin C helps to lower the non-specific binding of troponin I to the filter membranes.

TABLE 10

Filter	Fraction of troponin I recovered	
	Without troponin C	With troponin C
CytoSep	0.03 ± 0.03	0.15 ± 0.04
Glass Fiber	0.00 ± 0.03	0.09 ± 0.04

EXAMPLE 19

Use of Proteins of High Isoelectric Point to Prevent Non-specific Binding of Troponin I to Filter Membranes

A blood filter (CytoSep 1.5 cm×3.0 cm) was soaked for 16 hours at room temperature in solutions of deionized water containing 1 mg/ml of the proteins listed in Table 11. The filters were rinsed once with deionized water and dried for 2 hours at 35° C. Oxidized human cardiac troponin I (Bio-tech International, Inc.) at 100 ng/ml in human serum (Hybritech, Inc., San Diego, Calif.) also containing 2 mM added calcium chloride was passed through the filters as described in Example 18. The amount of troponin I recovered from the filters was determined by assay (Protocol B, Example 2) using TRI-7F81 anti-troponin I and biotinylated anti troponin I peptide 4 specific antibodies with added binding inhibitors (Example 11). The data in Table 11 are expressed as the fraction of troponin I recovered, which was determined as described in Example 18.

The results show that Melittin substantially reduces the non-specific binding of troponin I to the blood filter, whereas the other proteins had little effect at the concentrations tested.

TABLE 11

Protein	Fraction of troponin I recovered
No addition	0.16
Cytochrome C	0.11
Casein	0.16
Melittin	0.72
Lysozyme	0.15
Non-fat dried milk	0.08

EXAMPLE 20

Immunoassay of Ternary Troponin Complex using TRI-7 F81 anti Troponin I and Biotinylated anti Troponin I Peptide 4 Antibodies

The purified ternary troponin complex (Bio-tech International, Inc.) was assayed for troponin I as described in Example 13, except that the title antibody pair was used in the immunoassay and the aliquot of the reaction mixture of the antibodies with the troponin sample was taken three hours after the antibodies were added to the troponin. The fractional assay response was 0.16 in the absence of binding inhibitors and 0.49 in the presence of binding inhibitors.

The data show that the title antibody pair recognizes troponin I in the ternary complex poorly. In example 10, it was shown that the presence of troponin C without Troponin T had little effect on the recognition of the title antibody pair

for troponin I. Thus, the title antibody pair can bind to troponin I present in the binary complex with troponin C better than it can bind to Troponin I present in the ternary complex.

EXAMPLE 21

Immunoassay of Troponin I in Plasma From a Patient With Confirmed Myocardial Infarction Using TRI-7 F81 Anti Troponin I Conjugate and Biotinylated Anti Troponin I Peptide 1 Specific Antibodies

Frozen plasma from a patient with a confirmed myocardial infarction was thawed and diluted in human serum (Hybritech Inc., San Diego, Calif.) also containing 0.5M added sodium chloride and assayed for troponin I with the title antibody pair using Protocol A, Example 2. Oxidized purified human cardiac troponin I (University of Birmingham) was assayed and used to construct a standard curve from which the troponin I in the plasma was determined. The neat plasma sample was refrozen in a -70° C. freezer after being on ice for several hours. The frozen plasma was rethawed at room temperature and was assayed using the same protocol as described above except the plasma and troponin I standards were diluted into assay buffer. The standards exhibited almost identical assay slopes when diluted in serum (first assay) or assay buffer (second assay). The neat plasma sample was further incubated at 4° C. for the times indicated in Table 12 and reassayed in assay buffer.

The data show a substantial increase in the assayed troponin I concentration after a freeze/thaw cycle and after incubation at 4° C. This assay instability may be associated with the opening up or partial unraveling of the ternary troponin complex by the freeze/thaw cycle.

TABLE 12

Time of assay	Assayed concentration of troponin I
After first thaw of plasma	284 ng/ml
2 hours after second thaw of plasma	>800 ng/ml
19 hours after second thaw of plasma	1760 ng/ml
90 hours after second thaw of plasma	2300 ng/ml

We claim:

1. An assay for determining the presence or amount of a cardiac specific free and complexed isoform of troponin, said complex comprising isoforms selected from the group consisting of troponin I, troponin C and troponin T, comprising the steps of performing an immunoassay with antibody specific to a cardiac troponin which binds a complex comprising two or more of said isoform and another of said isoform of troponin in said sample, wherein the results of said immunoassay indicate the total amount of said free and complexed cardiac specific troponin isoform in said blood sample.

2. The method of claim 1, wherein the assay comprises a sandwich immunoassay.

3. The method of claim 1 wherein the performing step comprises providing:

- a conjugate antibody which binds to a cardiac specific isoform of troponin where the conjugate antibody is conjugated to a signal development element; or,

provides a solid phase antibody which binds to a cardiac specific isoform of troponin.

4. The method of claim 1 wherein the performing step comprises providing:

- 5 a conjugate antibody which binds to a cardiac specific isoform of troponin where the conjugate antibody is conjugated to a signal development element; and, a solid phase antibody which binds to a cardiac specific isoform of troponin.

5. The method of claim 1 wherein the performing step comprises providing an antibody which binds:

- complexes comprising the cardiac specific isoform of troponin I oxidized; complexes comprising the cardiac specific isoform of troponin I reduced, or, complexes comprising the cardiac specific isoform of troponin I oxidized, and complexes comprising the cardiac specific isoform of troponin I reduced.

6. An immunoassay method for determining the presence or amount of complexes of a cardiac specific isoform of troponin in a blood sample, said method comprising performing step a) or b), and steps c) through e) as follows:

- a) providing an antibody conjugate comprising an antibody which specifically binds to the cardiac specific isoform of troponin in a complex comprising said isoform and another isoform of troponin, and a means for generating a signal; or, b) providing an antibody conjugate comprising an antibody which specifically binds to the cardiac specific isoform of troponin in a binary complex comprising said isoform and another isoform of troponin, and a means for generating a signal; and, providing an antibody conjugate comprising an antibody which specifically binds to a cardiac specific isoform of troponin I or T, or troponin C in a ternary complex of three isoforms of troponin, and a means for generating a signal; and, c) contacting said sample with said antibody conjugate of step a), or with the antibody conjugates of step b),

whereby the antibody conjugates bind to the troponin complexes present in the sample;

- d) contacting the sample with a solid phase to which is bound at least one capture antibody which specifically binds said troponin complexes which are bound to antibody conjugate or bound to troponin complexes which is bound to antibody conjugate following binding to the capture antibody,

whereby a detectable signal is present from the antibody conjugate upon binding to the solid phase of cardiac specific troponin complexes which are or become bound to antibody conjugate; and,

- e) relating the detectable signal to the total amount of the cardiac specific isoform of troponin in said sample.

7. An immunoassay method for determining the presence or amount of complexes comprising two or more isoforms of troponin, of any cardiac specific isoform of troponin in a blood sample, said method comprising performing step a) or b), and steps c) through f) as follows:

- a) providing an antibody conjugate comprising an antibody which specifically binds to troponin complexes comprising two or more troponin isoforms, and a means for generating a signal; or,

- b) providing an antibody conjugate comprising an antibody which specifically binds to troponin in binary complexes comprising two troponin isoforms, and a means for generating a signal; or, providing an anti-

body conjugate comprising an antibody which specifically binds to ternary complexes comprising three troponin isoforms, and a means for generating a signal; and,

c) contacting said sample with said antibody conjugate of step a), or with the antibody conjugate(s) of step b), whereby the antibody conjugates bind to troponin complexes present in the sample;

d) contacting the sample with a solid phase to which is bound at least one capture antibody which specifically binds a cardiac specific isoform of troponin in complexes which are bound to antibody conjugate, or,

with a solid phase to which is bound at least one capture antibody which specifically binds to a cardiac specific isoform of troponin which is bound to antibody conjugate after binding to capture antibody,

whereby a detectable signal is present from the antibody conjugate upon binding to the solid phase of the cardiac specific troponin isoform which is or becomes bound to antibody conjugate; and,

e) relating the detectable signal to the presence or amount of the cardiac specific isoform of troponin in said sample.

8. An assay for determining the presence or amount of all free cardiac specific troponin isoforms and complexed isoforms of cardiac troponin in a blood sample, said complexes comprising isoforms selected from the group consisting of troponin I, troponin C and troponin T, comprising the steps of: performing an immunoassay with antibody specific to cardiac troponin which binds a complex comprising two or more of said isoforms of troponin in said sample, wherein the results of said immunoassay indicate the total amount of said free and complexed cardiac specific troponin isoforms in said blood sample.

9. An assay for determining the presence or amount of complexed isoforms of cardiac troponin in a blood sample, said complexes comprising isoforms selected from the group consisting of troponin I, troponin C and troponin T, comprising the steps of: performing an immunoassay with antibody specific to cardiac troponin which binds a complex comprising two or more of said isoforms of troponin in said sample wherein the results of said immunoassay indicate the total amount of said complexed cardiac specific isoforms in said blood sample.

10. The assay of claim 8 or 9, wherein said immunoassay provides a quantitative measurement for each of said isoform.

11. The assay of claim 1 wherein said free and complexed cardiac specific isoform is troponin I.

12. The assay of claim 1 wherein said free and complexed cardiac specific isoform is troponin T.

13. The assay of claim 1 wherein said assay provides a method for diagnosing a myocardial infarction, and said method includes correlating the results of said immunoassay with the existence of complexed troponin in said sample.

14. The assay of claim 1, 8 or 9, wherein said antibody comprises a plurality of monoclonal antibodies.

15. An assay for determining the presence or amount of a free and complexed cardiac specific isoform of troponin in a blood sample, said complex comprising isoforms selected from the group consisting of troponin I, troponin C and troponin T, comprising the steps of performing an immunoassay with antibody specific to cardiac troponin which binds said isoform and a complex comprising two or more of said isoform and another of said isoform of troponin in said sample, wherein the results of said immunoassay indicate the total amount of said free and complexed cardiac specific troponin isoform in said blood sample.

16. An assay for determining the presence or amount of a free and complexed cardiac specific isoform of troponin in a blood sample, said complex comprising isoforms selected from the group consisting of troponin I, troponin C and troponin T, comprising the steps of performing an immunoassay with antibody specific to cardiac troponin I and/or T which binds a complex comprising two or more of said isoform and another of said isoform of troponin in said sample, wherein the results of said immunoassay indicate the total amount of said free and complexed cardiac specific troponin isoform in said blood sample.

17. An assay for determining the presence or amount of a free and complexed cardiac specific troponin I in a blood sample, said complex comprising isoforms selected from the group consisting of troponin I, troponin C and troponin T, comprising the steps of performing an immunoassay with antibody specific to cardiac troponin I which binds a complex comprising two or more isoforms, one of which is troponin I and another of said isoform of troponin in said sample, wherein the results of said immunoassay indicate the total amount of said free and complexed cardiac specific troponin I isoform in said blood sample.

18. An assay for determining the presence or amount of a free and complexed cardiac specific troponin T in a blood sample, said complex comprising isoforms selected from the group consisting of troponin I, troponin C and troponin T, comprising the steps of performing an immunoassay with antibody specific to cardiac troponin T which binds a complex comprising two or more isoforms, one of which is troponin T and another of said isoform of troponin in said sample, wherein the results of said immunoassay indicate the total amount of said free and complexed cardiac specific troponin T isoform in said blood sample.

* * * * *